# Shop RANN USA Product Catalog 2025



Nothing beats the comfort and reliability of the Rainin LTS® LiteTouch® System

# **NEW Products to Enhance Your Workflows**



#### NEW

## More Depth, Breadth with MicroPro

Our new 1250  $\mu$ L model greatly expands MicroPro's volume range, covering practically any 96-, 384- as well as deep-well plate protocol | See page 24

mt.com/SR-MP

#### NEW

#### MicroPro 300 µL Tips in Strips

Ideal for serial dilution and/or partial plate filling. Racked tips are available for all models: 20, 300 and 1250 I See pages 25 and 30



# **NanoRep Repeater**

Feature packed. No Touch Off<sup>™</sup> dispense eliminates a source of cross contamination. 100 nL - 50 mL | See page 18

Pipettes	
LTS® LiteTouch® System	4
Single channel manual and electronic XLS+™	6
Multichannel manual and electronic XLS+	10
Adjustable Spacer XLS®	12
Rainin Classic™	14
XLS+ Starter Kits	15
Pipette performance specifications	16
NanoRep™ electronic repeater pipette	18
Rainin Pos-D™ positive displacement pipette	20
Positive vs. air displacement pipetting	21
High Throughput Pipetting	
96-Channel pipetting	22
MicroPro™ portable electronic system	24
BenchSmart™ 96 semi-automated system	26
Liquidator® 96 manual system	28
Accessories	29
Tips for 96-channel platforms	30
Specifications for 96-channel platforms	31
Liquid Handling – Specialty Products	
QuickFlow liquid aspiration	32
Rainin SP+ serological pipette controller	33
Rainin serological pipettes	33
AutoRep™ S manual repeater pipette	34
Disp-X <sup>™</sup> bottle top dispenser	35
Pipette Service	
Pipette service and service plans	36
Calibration vouchers	42
At Your Service! We Are Here to Help	43
GPP Good Pipetting Practice	44
Pipette Tips	
Reduce plastic waste with Rainin	46
Rainin EarthRack	48
TerraRack® tips	50
Green-Pak® SpaceSaver® and Green-Pak tips	52
StableStak™ StableRak™ and bulk tips	54
Racked tips	56
Low-retention tips	58
Specialty tips	60
BioClean Ultra™	62
Pipette Asset Management	64
SmartStand	65
SmartCheck™	66
PipetteX™ asset management software	67


# The Rainin LTS LiteTouch System

# Because Pipetting Shouldn't Hurt

Rainin revolutionized pipetting with the LiteTouch® System. LTS® not only virtually eliminates the harmful forces associated with tip loading and ejection, but also improves experimental reproducibility by eliminating inaccuracies due to imperfect seals.

#### Why switch to LTS?


- · Reduce your risk of pipetting-related RSI
- · Effortless tip insertion and ejection
- Perfect seals across all channels
- Pipette all day, every day without pain or fatigue

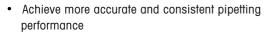



Meet LTS in 90 seconds! mt.com/LTS-video

#### Good ergonomics = better, more consistent results

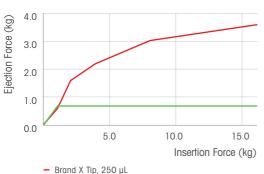


Rainin LTS Tip, 250 µL

Rainin LTS Tip, 250 µL

Accuracy requires consistency, yet consistency is difficult to maintain when a user's hands are fatigued or in pain. Pipetting is highly repetitive and the effects of pipetting forces – jamming tips onto shafts, then struggling to eject them – are cumulative, hence the documented relationship between pipetting accuracy and the effects of pipetting forces.

At 0.6 kg, LTS all but eliminates tip loading/ejecting forces, so users face none of the debilitating issues one experiences when working with standard tips.



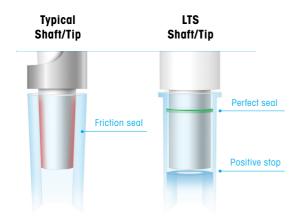

#### Effortless tip loading and ejection

Achieving a reliable seal with standard tips relies heavily on the strength (and persistence) of the user. Conventional pipette shaft-ends and tips are cone-shaped, so a tight seal relies on the strength of the user and pliability of the tip material. What's more, tip geometries vary from manufacturer to manufacturer, affecting the seal consistency and the accuracy of dispensed volumes.

LTS shafts and tips are cylindrical with a unique sealing ring inside the tip. Tips slide easily onto the shaft, stop when properly seated and make a perfect seal every time. Once seated, LTS tips won't fall off, yet slide off easily when pushed by the ejector.

#### Are you at risk?

- Pipetting is one of the most repetitive tasks in the laboratory.
- More than 30% of all workplace injuries are musculoskeletal disorders, including repetitive strain injuries and cumulative trauma disorders (injuries associated with repetition, excessive force and poor posture).
- Pipettes without LTS require up to 12 kg\* of force to complete a pipetting cycle.
- Forces should not exceed 30% of maximum average strength capacity:



- You are at risk if:
  - You pipette more than an hour a day or 300 hours a year.
  - You use non-LTS pipettes that can require up to 12 kg\* of force.
- Prevention is the key to reducing musculoskeletal disorders and repetitive strain injuries.

#### \*What is 12 kg of force equal to?

Driving a stake into the ground with a sledgehammer

Turning a manual can opener to open a can of food Picking up a small child

Download our pipetting ergonomics white paper mt.com/pipetting-ergo

# Rainin Pipet-Lite XLS+

# Performance You Can Feel

Rainin Pipet-Lite® XLS+ raises the bar on comfort, performance and control. Many volume sizes available in both LTS us and universal wy formats.

#### Smooth operation

The sealing system on XLS+ pipettes delivers a plunger stroke that is smooth, precise and requires less force than before. Perfect for pipetting applications that require delicate control!

#### Highly accurate

Built to last and manufactured with premium quality materials, XLS+ pipettes are designed to deliver superior performance, year after year.

#### **Autoclavable**

XLS+ single-channel pipettes are designed for easy autoclaving where it counts: the complete liquid end, including the seals and piston. Autoclaving won't compromise performance or feel.

#### Total comfort •

With sure-fit handles, contoured fingerhooks, lighter springs and low-drag seals, you can comfortably use Pipet-Lite XLS+ pipettes all day.



The tip ejector is made from a high density polymer that is corrosion free, autoclavable and highly resistant to chemicals. The ejector slides off in a single easy motion for cleaning.



## Need guaranteed volumes?

With fixed volume Pipet-Lite XLS+ you can set your pipette to your exact needs. Available in many volume sizes, LTS and UNV.

Contact your sales representative: 800 4 RAININ (800-472-4646) or is@rainin.com

ĺ	HNV
ı	OIAA

	CATALOG No.	VOLUME RANGE	PRICE
	17014413	0.1–2 μL	\$ 598.00
	17014409	0.5–10 μL	598.00
	17014412	2-20 µL	598.00
ts1	17014408	10-100 μL	598.00
Pipet-Lite XLS+	17014411	20–200 μL	598.00
et-Li	17014414	20-300 μL	598.00
Pip	17014407	100–1000 μL	598.00
	17014410	200–2000 μL	598.00
	17011801	500-5000 μL	598.00
	17011795	1-10 mL	598.00

Universal pipettes use universal tips. For pipette specifications, see page 16.

#### Advanced asset management

drift while allowing for rapid volume adjustment, even while wearing gloves.

The RFID tag on board every Rainin XLS+ pipette stores calibration data that is electronically read by the Rainin SmartStand and PipetteX software. See page 67 for details.



ITC

LIS			
	CATALOG No.	VOLUME RANGE	PRICE
	17014393	0.1–2 μL	\$ 598.00
	17014388	0.5–10 μL	598.00
	17014392	2-20 µL	598.00
+	17014384	10-100 μL	598.00
XLS	17014391	20-200 μL	598.00
Life	17014405	20-300 μL	598.00
Pipet-Lite XLS+	17014382	100-1000 μL	598.00
•	17014390	200-2000 μL	598.00
	17011790	500-5000 μL	598.00
	17011783	1-10 mL	598.00
	17011788	2-20 mL	598.00

LTS pipettes use LTS tips.
For pipette specifications, see page 16.

## You Deserve LTS

Happy Hands Yield Better Results!

Minimize fatigue and risk of repetitive strain injury (RSI) with the Rainin LTS® LiteTouch® tip ejection system.



The Rainin LiteTouch System dramatically reduces the amount of force required to load and eject tips. LTS pipettes work in concert with LTS tips to deliver a perfect seal every time.

- LTS tips slide easily over the shaft to create a firm, airtight fit with minimal force.
- A thin sealing ring near the top of an LTS tip virtually eliminates friction between the tip and shaft, minimizing ejection force.
- Positive stop prevents over-insertion of the tip and ensures an airtight seal.

#### mt.com/SR-LTS

# Rainin E4 XLS+

# Simple, Versatile and Secure

Our easiest to use, most richly featured and secure electronic pipette. Many volume sizes available in both LTS us and universal work formats.

#### The E4 XLS+ is easily configured to be as simple or elaborate as your work requires.

- The large color screen, joystick control and graphic interface make moving between functions and operating the pipette easy and effortless.
- Highly customizable functions with multiple protocols that can be stored onboard for later use.

Your data is everything – trust the E4 XLS+ to deliver!



#### • A "smart" charging stand!

SmartStand charges up to four E4s simultaneously and offers a simple yet sophisticated way to track and manage every XLS brand pipette in your lab. See page 65.

#### **Performance**

An advanced microprocessor-guided stepper motor has 4,000 discrete steps for extremely accurate liquid dispensing.

#### The app master

For serial dilutions, filling plates, pipetting unusual liquids and performing complex pipetting steps, the E4's versatile range of pipetting modes offers maximum convenience.

#### Save your protocols

Store multiple custom protocols in the onboard memory. Simultaneously review and select your protocols with a simple movement of the joystick. No more re-entering protocols over and over again!

#### Highly configurable

With E4's Admin Mode, you can configure the pipette any way you want. For convenience or security, remove unwanted modes or lock in preferred settings.





Minimize deviations from your pipetting SOPs by locking protocol settings. Onboard service data is extremely secure and quick to access. Avoid calibration compliance errors by password protecting access to service interval alarms.

mt.com/SR-E4

TC	
113	

LIO			
	CATALOG No.	VOLUME RANGE	PRICE
	17014484	0.5–10 μL	\$ 1,346.00
	17014487	2-20 μL	1,346.00
	17014483	10-100 μL	1,346.00
	17014486	20-200 μL	1,346.00
+ST	17014488	20-300 μL	1,346.00
E4 XLS+	17014482	100–1000 μL	1,346.00
	17014485	200-2000 μL	1,346.00
	17012312	500–5000 μL	1,346.00
	17012313	1-10 mL	1,346.00
	17012314	2-20 mL	1,346.00

LTS pipettes use LTS tips. Li-Ion battery and wall power supply are included. For more information about LTS, see page 4. For pipette specifications, see page 16.

	CATALOG No.	VOLUME RANGE	PRICE
	17014491	0.5-10 μL	\$ 1,346.00
	17014494	2-20 μL	1,346.00
	17014490	10–100 μL	1,346.00
<b>.</b>	17014493	20-200 μL	1,346.00
E4 XLS+	17014495	20–300 μL	1,346.00
7	17014489	100–1000 μL	1,346.00
	17014492	200–2000 μL	1,346.00
	17012353	500–5000 μL	1,346.00
	17012354	1-10 mL	1,346.00

Universal pipettes use universal tips. Li-Ion battery and wall power supply are included. For pipette specifications, see page 16.

# **Rainin XLS+ Multichannels**

# Performance You Can Feel!

8 or 12 channels; volumes from 0.5–1200  $\mu L.$  Available only in LTS.  $_{\mbox{\scriptsize LTS}}$ 

#### **E4 XLS+ Electronic Multichannel**

# Nothing compares with the functionality, speed and accuracy of E4 XLS+

- Same functionality and ease of use as single-channel models.
- E4 does the work with Multidispense and Autopace.
- Consistent sample loading across all channels with low tip ejection forces.


LTS			
	CATALOG No.	VOLUME RANGE	PRICE
	8-Channel		
	17013792	0.5–10 μL	\$ 2,340.00
	17013793	2-20 µL	2,340.00
nnel	17013794	5–50 μL	2,340.00
icha	17013795	20-200 μL	2,340.00
Mult	17013796	20–300 μL	2,340.00
nic	17014498	100–1200 μL	2,458.00
ectro	12-Channel		
Ä	17013797	0.5-10 μL	\$ 2,974.00
хГS-	17013798	2-20 µL	2,678.00
E4	17013799	5–50 μL	2,678.00
	17013800	20-200 μL	2,678.00
	17013801	20-300 μL	2,678.00
	17014499	100–1200 μL	3,596.00
E4 XLS+ Electronic Multichanne	12-Channel 17013797 17013798 17013799 17013800 17013801	0.5–10 μL 2–20 μL 5–50 μL 20–200 μL 20–300 μL	\$ 2,974.00 2,678.00 2,678.00 2,678.00 2,678.00

LTS pipettes use LTS tips. Li-lon battery and wall power supply are included. For more information about LTS, see page 4. For pipette specifications, see page 17.




# Industry-leading Channel-to-Channel Consistency

There's a reason Rainin multichannel pipettes are the market leader: ease-of-use and exceptional channel-to-channel consistency.

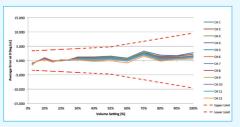


Rainin XLS+ 1200 µL



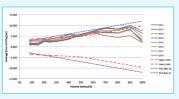


Competitor A • 1200 µL

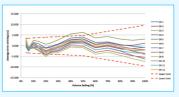

Competitor B • 1200 µL

VS.

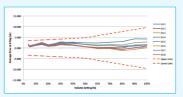

Competitor C • 1200 µL

# Pipet-Lite XLS+ Manual Multichannel Ergonomic design and exceptional performance

- Balanced, lightweight, reduced loading and ejection forces.
- Easy quick-lock volume setting no drift
- · Rugged and durable.

LTS	CATALOG No.	VOLUME RANGE	PRICE
	8-Channel		
	17013802	0.5-10 μL	\$ 1,586.00
<u>e</u>	17013803	2-20 μL	1,367.00
han	17013804	5-50 μL	1,367.00
I I I	17013805	20-200 μL	1,367.00
Ē	17013806	20–300 μL	1,367.00
XLS+ Manual Multichanne	17014496	100-1200 μL	1,806.00
<b>≥</b>	12-Channel		
	17013807	0.5-10 μL	\$ 1,816.00
Pipet-Lite	17013808	2-20 μL	1,708.00
ipet	17013809	5-50 μL	1,708.00
Δ.	17013810	20–200 μL	1,708.00
	17013811	20–300 μL	1,708.00
	17014497	100-1200 μL	2,566.00

LTS pipettes use LTS tips.

For more information about LTS, see page 4.

For pipette specifications, see page 17.

# **Rainin Adjustable Spacer Pipettes**

# Tubes to Plates with Just a Twist!

Adjust from standard 9 mm spacing to 14 mm or 19 mm; manual or electronic 6- and 8-channel versions. Available only in LTS. [115]

### Pipet-Lite XLS Adjustable Manual Multichannel

One-hand adjustable spacing.

#### Fast, convenient

Quickly change from format to format in a single step

#### Nozzle spacing ranges

6-channel: 9-19 mm 8-channel: 9-14 mm

#### Multiple formats

Use with 24-, 48- and 96-well plates and a variety of tube formats

#### **Applications**

Ideal for genomic, proteomic and cellular assay applications

LTS	CATALOG No.	VOLUME RANGE	PRICE
ole .	6-Channel		
ıstal	17011841	20–300 μL	\$ 2,264.00
Adjı	17011840	100-1200 μL	2,264.00
Pipet-Lite XLS Adjustable	8-Channel		
Lite	17011844	5–50 μL	\$ 2,264.00
pet-	17011843	20–300 μL	2,264.00
<u>~</u>	17011842	100-1200 μL	2,264.00

LTS pipettes use LTS tips. For more information about LTS, see page 4. For pipette specifications, see page 17.



## **E4 XLS Adjustable Electronic Multichannel**

Different formats? No problem!

#### Easy, accurate and repeatable spacing

Just set the limiter dial and quickly move between formats by twisting the spacing control. Marks on the liquid end give fast visual confirmation of nozzle spacing. Ideal for routine micro-well work when you need to change formats quickly.

LTS	CATALOG No.	VOLUME RANGE	PRICE
	6-Channel		
ple	17012327	20–300 μL	\$ 3,830.00
E4 XLS Adjustable	17012328	100-1200 μL	3,830.00
Adj	8-Channel		
XLS	17012329	5–50 μL	\$ 3,826.00
E4	17012330	20-300 μL	3,830.00
	17012331	100–1200 μL	3,826.00

LTS pipettes use LTS tips. Li-lon battery and wall power supply are included. For more information about LTS, see page 4. For pipette specifications, see page 17.

# Reservoirs for Multichannel Pipettes

- · Capacity 75 mL
- · Available with and without lid
- · Great addition for repeat dispensing

CATALOG No.		QUANTITY	PRICE	
	30852500	Reservoir with lid, non-sterilized	25/pack	\$ 130.80
	30852501	Reservoir without lid, non-sterilized	50/pack	151.80





# **Rainin Classic**

# Timeless Design, Exceptional Performance

An economical pipette with advanced ergonomic features not found on other pipettes in its class. Available in universal w formats.

An economical pipette with decades of expertise built in, Rainin Classic™ will deliver years of reproducible performance. Lighter spring forces, low-drag seals and a shock absorber make it easy on your hands.

#### **Durable**

Stainless steel and high-quality PVDF materials ensure durability, fewer out-of-tolerance failures and low cost of ownership.

#### Lower plunger force

The Rainin Classic's plunger force is 25% less than similar pipettes.

#### Comfort

The fingerhook allows a lighter grip during pipetting cycles, for improved ergonomics.

#### Great value – lower prices!

UNV			
	CATALOG No.	VOLUME RANGE	PRICE
	17008648	0.1-2 µL	\$ 415.00
	17008649	0.5–10 μL	415.00
sic	17008650	2-20 μL	415.00
Clas	17008651	10-100 μL	415.00
Rainin Classic	17008652	20-200 μL	415.00
Rai	17008653	100–1000 μL	415.00
	17008654	500–5000 μL	415.00
	17008655	1-10 mL	415.00
Starter Kit	17008708	Starter kit: PR-20, PR-200, PR-1000 and accessories	\$ 1,042.00





# **Pipette Starter Kits**

# Upgrade Your Lab with New Rainin XLS+ Pipettes!

Available with SmartStand or Hang-Ups

Starter Kits with SmartStand include PipetteX software for easy pipette management

 Manage calibration and service for all of your pipettes

Single Channel

- View pipette calibration status right at your bench
- Charge your E4 XLS+ electronic pipettes



Multichannel



Single Channel

				Manual				Manual		E	Electronic						
		/			/ ,	/ ,	/ ,		200 Schonney	20	12-channe	20	,	/ ,	//	Hano Solvoie	8
LTS	UNV		, / .	s / 3	20°	¥ / ¿	300°	30 M	% %	¥/:	× / ,	¥/:	20°	¥ / §	Smorts		5
	CATALOG No.	\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\		20 M	/ 00	/ 0	, \ 0,0	100	000	20 K	0,00	20 M	200		/50	Hono	PRICE
	30386738	✓	✓	✓	✓										✓		\$ 2,040.00
	30579469	<b>✓</b>	✓	✓	<b>✓</b>											✓	1,760.00
	30386597	<b>✓</b>		✓	<b>✓</b>	✓									<b>√</b>		1,892.00
	30579367	<b>✓</b>		✓	<b>✓</b>	<b>✓</b>										✓	1,760.00
	30386599			✓	<b>✓</b>	✓									<b>√</b>		1,754.00
	30386737				<b>✓</b>	<b>✓</b>	✓								<b>√</b>		2,040.00
	30386734			✓	<b>✓</b>	<b>✓</b>		<b>√</b>							<b>√</b>		2,622.00
113	30386733			✓	<b>✓</b>	<b>✓</b>			<b>✓</b>						<b>√</b>		2,622.00
	30386736			✓	<b>✓</b>	<b>✓</b>				<b>√</b>					<b>√</b>		2,912.00
	30386735			✓	<b>✓</b>	<b>✓</b>					✓				<b>√</b>		2,912.00
	30386739							✓	✓						✓		2,326.00
	30456871			✓	<b>✓</b>	<b>✓</b>										✓	1,382.00
	30579472							<b>√</b>	✓							✓	2,162.00
	30386740									<b>✓</b>	<b>✓</b>				<b>√</b>		2,326.00
	30579471									<b>√</b>	✓					✓	2,162.00
	30386731											✓	✓	<b>√</b>	✓		3,928.00
	30386598	<b>✓</b>		<b>√</b>	<b>✓</b>	<b>✓</b>									<b>√</b>		\$ 1,892.00
	30386730			1	1	<b>✓</b>									<b>✓</b>		1,754.00
N N	30456872			<b>✓</b>	<b>✓</b>	<b>✓</b>										<b>✓</b>	1,382.00
	30386732											<b>✓</b>	1	<b>✓</b>	<b>✓</b>		3,928.00
	17008708*			✓	<b>✓</b>	<b>✓</b>										<b>✓</b>	1,042.00
	* Rainin Classi	С													,		

# **Pipette Performance Specifications**

# Specifications at 10%, 50% and 100% volume

# All Rainin single-channel models, including:

#### Manual pipettes

- Pipet-Lite XLS/XLS+
- Rainin Classic

#### **Electronic pipettes**

• E4 XLS/XLS+



		VO	L	INCREMENT µL		SYSTEMATI	C ERROR	RANDOM EI	ROR
					ω /				
		ME	1						
	PIPETTE VOLUME	Jun 101	MANULA	ELECTRO	/ %	(±) μL	/ %	(≤) μL	
	2 μL	0.2 μL	0.002	n/a	12	0.024	6.0	0.012	
		1 μL			2.7	0.027	1.3	0.013	
		2 μL			1.5	0.030	0.70	0.014	
	10 μL	1.0 µL	0.02	0.01	2.5	0.025	1.2	0.012	
		5.0 µL			1.5	0.075	0.60	0.030	
		10.0 µL			1.0	0.10	0.40	0.040	
	20 µL	2 μL	0.02	0.02	7.5	0.15	2.0	0.040	
		10 µL			1.5	0.15	0.50	0.050	
		20 µL			1.0	0.20	0.30	0.060	
	100 µL	10 µL	0.2	0.1	3.5	0.35	1.0	0.10	
		50 μL			0.80	0.40	0.24	0.12	
		100 µL			0.80	0.80	0.15	0.15	
	200 μL	20 µL	0.2	0.2	2.5	0.50	1.0	0.20	
tes		100 µL			0.80	0.80	0.25	0.25	
Single-Channel Pipettes		200 μL			0.80	1.6	0.15	0.30	
e P	300 µL	30 µL	0.5	0.2	2.5	0.75	1.0	0.30	
an		150 µL			0.80	1.2	0.25	0.375	
<u>ရ</u>		300 µL			0.80	2.4	0.15	0.45	
ngle	1000 µL	100 µL	2	1	3.0	3.0	0.60	0.60	
Si		500 μL			0.80	4.0	0.20	1.0	
		1000 µL			0.80	8.0	0.15	1.5	
	2000 μL	200 µL	2	2	3.0	6.0	0.60	1.2	
		1000 µL			0.80	8.0	0.20	2.0	
		2000 μL			0.80	16	0.12	2.4	
	5000 μL	500 μL	5	5	2.4	12	0.60	3.0	
		2500 µL			0.60	15	0.20	5.0	
		5000 μL			0.60	30	0.16	8.0	
	10 mL	1 mL	20	10	5.0	50	0.60	6.0	
		5 mL			1.0	50	0.20	10	
	00	10 mL	0.0	0.0	0.60	60	0.16	16	
	20 mL	2 mL	20	20	5.0	100	0.60	12	
		10 mL			1.0	100	0.20	20	
		20 mL			0.60	120	0.16	32	

# All Rainin multichannel models, including:

#### Manual pipettes

- Pipet-Lite XLS/XLS+
- Pipet-Lite XLS Adjustable Spacer

#### **Electronic pipettes**

- E4 XLS/XLS+
- E4 XLS Adjustable Spacer

		V	VOL   INCREMENT µL   SYSTEMATIC ERROR   MANUAL						RANDOM ERROR ELECTRONIC
	PIPETTE VOLUME	NOUME	MANUAL ELECTO	% %	(±) µL	%	(≤) µL	%	(≤) μL
	10 µL	1.0 µL	0.02/0.01	4.0	0.040	5.0	0.050	3.0	0.030
		5.0 µL		1.5	0.075	1.5	0.075	0.80	0.040
		10.0 μL		1.0	0.10	0.50	0.050	0.50	0.050
	20 µL	2 μL	0.02	7.5	0.15	4.0	0.080	3.0	0.060
		10 µL		1.5	0.15	1.5	0.15	1.0	0.10
		20 µL		1.0	0.20	0.50	0.10	0.30	0.060
ttes	50 μL	5 μL	0.05	3.5	0.175	2.0	0.10	1.5	0.075
ipe		25 µL		1.2	0.30	0.50	0.125	0.40	0.10
Multichannel Pipettes		50 μL		0.80	0.40	0.40	0.20	0.30	0.15
la la	200 µL	20 µL	0.2	2.5	0.50	1.2	0.24	1.0	0.20
= ic		100 µL		0.80	0.80	0.25	0.25	0.25	0.25
≥		200 μL		0.80	1.6	0.30	0.60	0.20	0.40
	300 µL	30 µL	0.5/0.2	2.5	0.75	1.5	0.45	1.0	0.30
		150 µL		0.80	1.2	0.30	0.45	0.25	0.375
		300 µL		0.80	2.4	0.25	0.75	0.20	0.60
	1200 µL	100 µL	2/1	3.6	3.6	1.0	1.0	0.60	0.60
		600 µL		0.80	4.8	0.30	1.8	0.20	1.2
		1200 µL		0.80	9.6	0.20	2.4	0.15	1.8





Specifications subject to change without notice.

# Rainin NanoRep

# Electronic Repeater Pipette

Automatically dispense up to 1,000 aliquots

- with no contact necessary from vessel to vessel.
- Lightweight, comfortable grip
- Transfer volatile, viscous and dense liquids
- No Touch Off<sup>™</sup> dispense
- Aliquots as low as 100 nL, with increments as low as 10 nL – down to 1/10,000 th the nominal volume of any syringe tip
- Up to 1,000 aliquots from one aspiration



	CATALOG No.	DESCRIPTION	PRICE
۵	30568171	Electronic NanoRep repeater pipette (includes: charging cable, Hang-Up and 10 mL syringe tip)	\$ 1,862.00
NanoRep	17012878	Wall power supply	81.00
P S	30584027	Hang-Ups (for pipettes with fingerhook)	95.65
	30575784	NanoRep stand	126.80



#### NanoRep syringe tips I Specifications

Tip Vol.	1%	Systema 10%	tic Error <b>50%</b>	100%
0.1	1 μL	10 μL	50 μL	100 μL
	± 8.0%	± 1.6%	± 0.90%	± 0.90%
	(0.080 μL)	(0.16 μL)	(0.45 μL)	(0.90 μL)
10	100 μL	1 mL	5 mL	10 mL
	± 2.0%	± 0.40%	± 0.40%	± 0.30%
	(2.0 μL)	(4.0 μL)	(20 μL)	(30 μL)
50	500 μL	5 mL	25 mL	50 mL
	± 2.0%	± 0.25%	± 0.25%	± 0.25%
	(10 μL)	(12.5 μL)	(62.5 μL)	(125 μL)

Tip Vol.	1%	Randor 10%	n Error <b>50%</b>	100%
0.1	1 μL	10 μL	50 μL	100 μL
	≤ 12%	≤ 2.5%	≤ 0.80%	≤ 0.50%
	(0.12 μL)	(0.25 μL)	(0.40 μL)	(0.50 μL)
10	100 μL	1 mL	5 mL	10 mL
	≤ 3.0%	≤ 0.35%	≤ 0.20%	≤ 0.15%
	(3.0 μL)	(3.5 µL)	(10 μL)	(15 μL)
50	500 μL	5 mL	25 mL	50 mL
	≤ 3.5%	≤ 0.45%	≤ 0.20%	≤ 0.10%
	(17.5 μL)	(22.5 μL)	(50 μL)	(50 μL)

Tip Volume

mL	Min. Aliquot* <b>µL</b>	Min.Incr.** µL
0.1	0.1	0.01
10	10	1
50	50	5

<sup>\*</sup>Min. Aliquot = Smallest dispensable volume/tip.

	CATALOG No.	MAX. VOLUME	QUANTITY		PRICE
Tips	30575705	0.1 mL	100 per pkg		\$ 183.80
	30575777	0.1 mL	25 per pkg	Sterilized	58.20
Syringe	30575707	10 mL	100 per pkg		183.80
	30575779	10 mL	25 per pkg	Sterilized	58.20
NanoRep	30575789	50 mL	25 per pkg		79.80
N P	30575781	50 mL	25 per pkg	Sterilized	91.60

mt.com/SR-NR

<sup>\*\*</sup>Min. Increment = Lowest incremental adjustment/tip.

# **Rainin Pos-D**

# Positive Displacement Pipetting

The right solution for liquids with high densities, viscosities or vapor pressures.

Rainin Pos- $D^{\text{TM}}$  positive displacement pipettes are ideal when pipetting liquids with high densities, viscosities or vapor pressures.

#### Relax!

The ergonomically designed Pos-D rests in your hand while you work.

#### Easy to load

Pre-assembled, sterilized disposable syringe tips are packed in convenient racks.

#### No cross contamination

Positive displacement prevents aerosol formation and protects samples from contamination.

		VOL	·	QT >		
	04741.00 N-	VOLUME RANGE	700	TIPS/CARTON	STED	PRICE
	17008575		/ 4		/ ~	\$ 571.00
Ś	17008575	0.5–10 μL				571.00
Pos-D Pipettes		3–25 μL				
P.	17008577	20–50 μL				571.00
Q-S	17008578	10–100 μL				571.00
2	17008579	50-250 μL				571.00
	17008580	100–1000 μL				571.00
	17008604	10 μL	3	180	<b>√</b>	\$ 116.20
	17008605	25 µL	3	180	✓	116.20
	17008606	50 μL	3	180	<b>✓</b>	116.20
SUO	17008607	100 μL	3	180	<b>√</b>	116.20
Pist	17008608	250 μL	3	180	<b>√</b>	116.20
Capillaries and Pistons	17008609	1000 μL	3	180	✓	116.20
ies	17012264	10 μL	3	180		105.40
ii ar	17012265	25 μL	3	180		105.40
Cap	17012266	50 μL	3	180		105.40
	17012267	100 μL	3	180		105.40
	17012268	250 μL	3	180		105.40
	17012135	1000 μL	3	180		105.40



# Why Consider Positive Displacement?


# A Quick Roundup

#### Ideal for wide range of liguids

Positive displacement outperforms air displacement in the handling of viscous and volatile liquids. Because they can be slow and difficult to move in and out of a pipette tip, viscous liquids can stretch and compress the air column inside an air displacement pipette, leading to measurement inaccuracies. Volatile liquids constantly emit vapor, so when aspirated into an air displacement pipette tip, air pressure builds up inside the pipette and can push liquid out of the tip without any dispense action initiated by the user. This can also lead to measurement inaccuracies.

#### Extreme accuracy

Positive displacement pipettes use syringe tips that contain an inner plastic piston extending the full length of the tip. At zero volume, the end of the inner piston is fully extended, closing off the opening of the tip. As a user aspirates liquid, the end of the piston rises, pulling liquid into the tip. With zero air between the piston and the liquid in a positive displacement system, the pressures that viscous and volatile liquids exert on air do not come into play, and extreme accuracy is preserved.

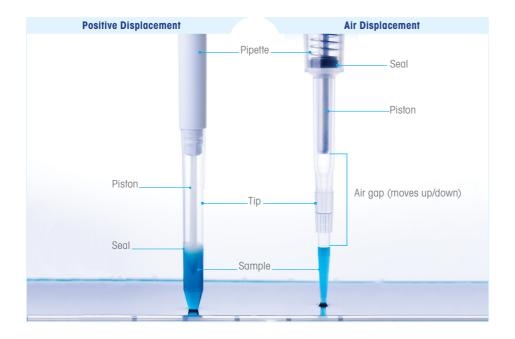

# 96-Channel Pipetting

# Faster, More Accurate Workflows

By making it possible to perform experiments more efficiently, accurately and reliably, 96-channel pipettors have become indispensable tools in many life science laboratories. They can significantly reduce the time and effort required to perform experiments, while providing higher accuracy and reproducibility. They are particularly useful in applications that involve large-scale screenings, such as drug discovery, ELISAs, proteomics and cell-based assays.


Unlike single- and multichannel pipettes, which are more dependent on user's technique, manual and semi-automated 96-channel pipettors are ideal for research and scale-up, and can often provide more reliable and reproducible results. They also dispense small volumes of liquid with high accuracy, which is important in many applications that require precise measurements.

Rainin offers 96-channel pipettors to meet your exact need.

- MicroPro our most compact, affordable and advanced semi-automated 96-channel pipettor.
- BenchSmart 96 with four plate-trays, ideal for labs that do a lot of multi-plate work.
- Liquidator 96 perfect for researchers who want to manually control every aspect of their workflow.


	MicroPro	BenchSmart 96	Liquidator 96
96/384-Well Plate Applications	•	•	•
Interchangeable Heads	•	•	•
Suitable for BSL/PCR Cabinets	•	•	•
Small Footprint	•	•	•
Low Cost	•	•	•
Volume Range	0.5-20 μL 5-300 μL 50-1250 μL	0.5-20 μL 5-200 μL 100-1000 μL	0.5-20 μL 5-200 μL

# How can a 96-channel pipettor speed your workflow?

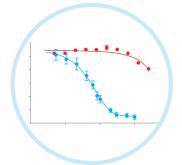


#### High-Throughput Screening

A 96-channel pipettor enables you to significantly increase your screening capacity and streamline workflows. You can quickly and accurately screen large numbers of samples for various applications, such as genomic, drug screening or proteomics research.

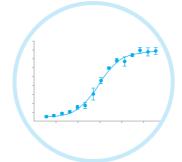


#### **ELISA Assays**

Enzyme-Linked Immunosorbent Assays (ELISA) are a common method for detecting and quantifying proteins or small molecules in biological samples. Using a 96-channel pipettor, you can simultaneously transfer the reagents to multiple wells or wash wells, significantly reducing the time and effort needed to perform an ELISA assay.

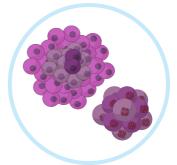


#### Cell-Based Assays

Cell-based assays are used to study the behavior of cells in response to various stimuli. Using a 96-channel pipettor, you can quickly and accurately transfer cells, reagents and other components to multiple wells, enabling you to perform high-throughput cell-based assays.

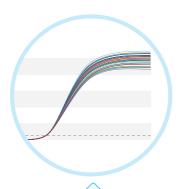

#### PCR/qPCR Amplification

Polymerase Chain Reaction (PCR) is a method used to amplify DNA and RNA sequences. Using a 96-channel pipettor, you can quickly and accurately transfer the reagents needed for PCR amplification, streamlining the process and increasing the throughput.

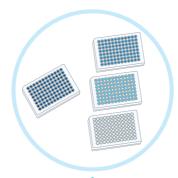

#### Plate Replication

Plate replication is commonly used in life science research to create multiple copies of a plate for analysis. Using a 96-channel pipettor, you can quickly and accurately transfer the contents of a mother plate to multiple daughter plates, increasing efficiency and reducing the chance of errors.



#### **DNA Sequencing**

DNA sequencing is a critical technique used in genomics research. Using a 96-channel pipettor, researchers can quickly and accurately accomplish the many sample preparation steps often involved prior to DNA sequencing.

## Rainin MicroPro

# Affordable, Portable 96-Channel Pipetting

MicroPro simplifies and streamlines 96- and 384-well microplate work, from filling plates to complex, multi-step workflows.

NEW 1250 µL model

Compact, highly affordable and easy to use, the MicroPro wirelessly connects to a high resolution touchscreen for an intuitive and user-friendly experience.

- 20 μL, 300 μL and 1250 μL models
- Weighs just around 13 pounds!
- Save/retrieve complex workflows
- · Intuitive touchscreen interface
- · 10x speed control





# Simplify your workflows

From the touchless tip ejection and programmable Pipetting Depth Recall™ to saving complex, multistep protocols, the MicroPro boosts productivity and improves reproducibility.



#### Your research assistant

MicroPro goes beyond simply adding and mixing reagents by offering convenient features such as aliquot dispensing, sample dilution with mixing and multi-aspirate with mixing.



#### Custom protocols

Custom protocols are easy to create and allow you to build and save a tailored sequence of instruction and liquid handling steps that exactly match the needs of your protocols.

mt.com/SR-MP



#### MicroPro Tips

The innovative and proprietary design of the MicroPro pipette tips and tip cartridges requires minimal force for loading all 96 tips. This significant force reduction allows MicroPro to be small and lightweight — currently the smallest benchtop on the market! MicroPro tips work with the Rainin MicroPro instrument only. They are available as low-retention tips in 20, 300 and 1250  $\mu L$  volume sizes, sterilized and sterilized with filter.



With our new MicroPro 300  $\mu$ L "Tips in Strips" you can perform serial dilutions (a series of dilutions related to each other by setting a constant dilution) and protocols involving partial plate filling (replicates within the same plate) with greater ease, speed and reproducibility.



#### MicroPro

	CATALOG No.	MAX VOL	DESCRIPTION	PRICE
	30835029	20 μL	0.5-20 µL MicroPro 20	\$ 13,470.00
MicroPro	30835028	300 µL	5-300 µL MicroPro 300	13,470.00
Micr	30932001	1250 µL	50-1250 μL MicroPro 1250	13,470.00
	30835030		Two-Tray Dual Nest, MPR 96/96	1,616.00

See page 30 for high-throughput MicroPro tips.

## Rainin BenchSmart 96

# Semi-Automated 96- and 384-Well Pipetting

The precision and reproducibility of automated liquid handling with the speed and flexibility of a manual system.

# The semi-automated BenchSmart™ 96 improves reproducibility by eliminating user variability when aspirating and dispensing.

- Semi-automated high-throughput 96- and 384-well plate processing for genomics, proteomics, and cellular-based assays.
- BenchSmart gives users complete and precise control over the position of the pipetting head with 96 channels.
- The four tray positions save time and reduce your risk of error by minimizing the need to swap out trays and reservoirs.





#### Smart high-throughput

With its large touchpad, BenchSmart makes it easy to design, save and retrieve any protocol, from simple one-step procedures to complex, multi-step experiments.



#### Three pipetting heads

Get the precision you demand with the 0.5-20  $\mu$ L head, the lab workhorse you need with the 5-200  $\mu$ L head, and the multi-dispensing capability of the 100-1000  $\mu$ L head.



#### Save time and money

BenchSmart 96 simplifies and speeds workflows by significantly reducing, if not eliminating, the need to swap out tips and reservoirs.

mt.com/SR-BST96




#### Speed your workflow

Compared to a two-tray system where reservoirs and trays need to be swapped out, a two solution protocol needs fewer liquid handling steps with BenchSmart.



4-tray BenchSmart 96

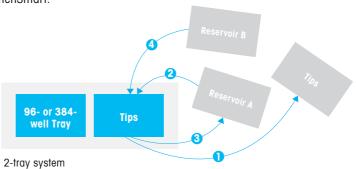

#### BenchSmart 96



	CATALOG No.	MAX VOL	DESCRIPTION	PRICE
	30296705	20 μL	0.5-20 µL BenchSmart 96, complete	\$ 26,940.00
	30296706	200 μL	5-200 µL BenchSmart 96, complete	26,940.00
96	30296707	1000 μL	100-1000 µL BenchSmart 96, complete	26,940.00
BST	30296708	20 μL	Pipetting head 0.5-20 $\mu$ L, BenchSmart 96	10,850.00
	30296709	200 μL	Pipetting head 5-200 µL, BenchSmart 96	10,850.00
	30296780	1000 μL	Pipetting head 100-1000 µL, BenchSmart 96	10,850.00
	17010394	-	384-well adapter plate stage, white	\$ 471.00
Accs.	17010791	_	384-well adapter plate stage, black	471.00
	30321288	-	Terminal Conversion Kit	529.00

See page 30 for high-throughput LTS tips. For reservoirs, plates, mats and strips, see page 29.

# **Rainin Liquidator 96**

# Manual High-Throughput Pipetting

Fully manual processing of 96- and 384-well plates for genomics, proteomics and cell-based assays.

Liquidator™ 96 streamlines application workflows for genomics, proteomics and cellular applications, simplifies plate replication, runs serial dilutions, performs plate washes or quickly moves samples between 96- and 384-well formats.

Count on Liquidator 96 to deliver the highest accuracy for PCR, qPCR, DNA sequencing and genotyping experiments. Generate precise 96-well ELISA, functional, cell viability or chemotaxis assays, all on one platform.


- Two volume ranges: 0.5-20 µL and 5-200 µL
- No power or programming: simply set the volume and pipette
- Four-position platform: process up to 3 plates at a time without stacking






#### Easy to use

Because Liquidator doesn't require programming or significant training, anyone can be an expert in no time.

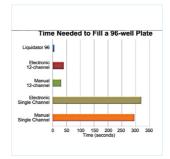

#### Save time

Liquidator processes 96-well plates in as few as six seconds each, without sacrificing data quality.

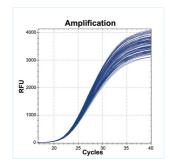

#### Reduce errors

By pipetting all 96 wells simultaneously, Liquidator eliminates the risk of skipping or repeating wells or rows.

mt.com/SR-LIQ96



# Liquidator 96



	CATALOG No.	MAX VOL	DESCRIPTION	PRICE
	17014207	20 µL	0.5-20 µL Benchtop Pipetting System	\$ 26,130.00
LIQ	17010335	200 μL	5-200 µL Benchtop Pipetting System	26,130.00

Liquidator 96 is supplied with height adjustment posts, a 384-well adapter plate stage and three racks of tips to get you started.

	17014270	-	Height adjustment posts, 0.5-20 µL Liquidator	\$ 88.40
Accs.	17010396	_	Height adjustment posts, 5-200 µL Liquidator	87.20
96 A	17010394	_	384-well adapter plate stage, white	471.00
FI 0	17010791	_	384-well adapter plate stage, black	471.00
	17011118	-	20 μL tip tray adapter, white	408.00

# Reservoirs, plates, mats and strips

Reusable SBS footprint polypropylene reservoirs

CATALUG	NO.	DESCRIPTION

PRICE

	17012602	Non-sterile	Low-profile 96 pyramidal bottoms	5-pack	\$ 75.00
	17012603	Sterile	Low-profile 96 pyramidal bottoms	5-pack, ind. wrap	108.80
	17012604	Non-sterile	Standard profile 96 pyramidal bottoms	5-pack	87.40
	17012605	Sterile	Standard profile 96 pyramidal bottoms	5-pack, ind. wrap	121.60
ဟ	17012608	Non-sterile	Low-profile 8-channel V-bottom	5-pack	75.00
Voir	17012609	Sterile	Low-profile 8-channel V-bottom	5-pack	108.80
Reservoirs	17012606	Non-sterile	Standard profile 8-channel V-bottom	5-pack	87.40
~	17012607	Sterile	Standard profile 8-channel V-bottom	5-pack, ind. wrap	121.60
	17012612	Non-sterile	Low-profile 12-channel V-bottom	5-pack	75.00
	17012613	Sterile	Low-profile 12-channel V-bottom	5-pack, ind. wrap	108.20
	17012610	Non-sterile	Standard profile 12-channel V-bottom	5-pack	86.80
	17012611	Sterile	Standard profile 12-channel V-bottom	5-pack, ind. wrap	121.60
	17012623	Non-sterile	2.2 mL 96-deepwell plate	5-pack	\$ 43.80
ries	17012624	Sterile	2.2 mL 96-deepwell plate	5-pack, ind. wrap	75.00
988	17012625	Non-sterile	Polymer sealing mat fits 96-deepwell plate	5-pack	43.80
Accessories	17012626	Sterile	Polymer sealing mat fits 96-deepwell plate	5-pack, ind. wrap	75.00
	17012630	Sterile	1.2 mL racked microtube strips (8x12)	box of 300	72.40

Deepwell plates, mats, tube strips, cap strips

# Tips for 96-Channel Platforms

High-Throughput	LTS	S Tips								
For use with BenchSmart 9			!	VOL		QT	Υ		FE	ATURES
pipettors only (1000 μL fo	r Be	nchSmart only).		Max VOLUME	. / 5	TIPS CARTON	STE	FILZ	1/ER 20/4.	PRICE
		CATALOG No.	DESCRIPTION	<u>z,</u>	67	/ 🗷	/ 5	/ 4		PRICE
		17011185	LQR LTS 20µL 960/10	20 μL	10	960				\$ 120.20
		17011186	LQR LTS 20µL S 960/10	20 μL	10	960	✓			139.60
		17011117	LQR LTS 20µL F 960/10	20 μL	10	960	✓	✓		200.40
		17010645	LQR LTS 200µL 960/10	200 µL	10	960				120.20
		17010647	LQR LTS 200µL S 960/10	200 μL	10	960	✓			139.60
	w	17010646	LQR LTS 200µL F 960/10	200 μL	10	960	✓	✓		200.40
	Tips in Racks	30281704	LQR LTS 1000µL 960/10	1000 μL	10	960				147.20
	교 드	30296781	LQR LTS 1000µL S 960/10	1000 μL	10	960	✓			168.60
	ips	30296782	LQR LTS 1000µL F 960/10	1000 μL	10	960	✓	<b>√</b>		238.60
		17014399	LQR LTS 20µL LS 960/10	20 µL	10	960	<b>√</b>		<b>✓</b>	168.60
		17014400	LQR LTS 20µL FL 960/10	20 μL	10	960	✓	<b>√</b>	<b>✓</b>	246.80
		17014401	LQR LTS 200µL LS 960/10	200 μL	10	960	✓		<b>✓</b>	168.60
		17014402	LQR LTS 200µL FL 960/10	200 μL	10	960	<b>√</b>	<b>√</b>	<b>✓</b>	246.80
		30296783	LQR LTS 1000µL LS 960/10	1000 μL	10	960	<b>√</b>		<b>✓</b>	168.60
		30296784	LQR LTS 1000µL FL 960/10	1000 μL	10	960	✓	✓	✓	246.80
		17011187	LQS LTS 20µL 960/10	20 μL	10	960				\$ 120.20
A STATE OF THE STA	s k k	17011287	LQS LTS 20µL S 960/10	20 µL	10	960	<b>✓</b>			139.60
	Stacks	17010648	LQS LTS 200µL 960/10	200 μL	10	960				120.20
	Lips in	17010649	LQS LTS 200µL S 960/10	200 μL	10	960	<b>√</b>			139.60
	Tips	17014403	LQS LTS 20µL LS 960/10	20 μL	10	960	<b>√</b>		<b>√</b>	168.60
		17014404	LQS LTS 200µL LS 960/10	200 μL	10	960	1		<b>✓</b>	168.60

#### MicroPro Tips in Racks



	63305398	MPR 20µL LS 4800/50	20 μL	50	4800	<b>√</b>		<b>√</b>	\$ 521.00
	63305397	MPR 20µL LF 4800/50	20 μL	50	4800	$\checkmark$	✓	<b>√</b>	718.00
so.	63305400	MPR 300µL LS 4800/50	300 µL	50	4800	$\checkmark$		<b>√</b>	484.00
Tips	63305399	MPR 300µL LF 4800/50	300 μL	50	4800	$\checkmark$	✓	<b>√</b>	698.00
MicroPro	63305538	MPR 1250µL LS 960/10	1250 µL	10	960	$\checkmark$		<b>√</b>	138.00
Aicr	63305537	MPR 1250µL LF 960/10	1250 µL	10	960	$\checkmark$	✓	✓	207.00
	Tips in Strips								
	63305551	MPR 300µL LS strips 960/10	300 µL	10	960	$\checkmark$		✓	\$ 120.60
	63305401	MPR 300µL LF strips 960/10	300 µL	10	960	$\checkmark$	✓	✓	238.00

# Specifications for 96-Channel Platforms

MODEL/MAX. VOLUME 20 μL

**Volume Increments** 

	MODEL/MAX. VOLUME	20 μL	200 μL	1000 μL
BenchSmart 96	HEAD SIZE	LOW VOLUME 0.5 – 20 μL	MID VOLUME 5 – 200 μL	HIGH VOLUME 100 – 1000 μL
	Systematic Error (per channel)	20 μL: ± 1.0% (0.20 μL) 10 μL: ± 1.2% (0.12 μL) 2 μL: ± 6.0% (0.12 μL) 1 μL: ± 12% (0.12 μL)	200 μL: ± 1.0% (2.0 μL) 100 μL: ± 1.0% (1.0 μL) 20 μL: ± 2.0% (0.40 μL) 5 μL: ± 5.0% (0.25 μL)	1000 μL: ± 1.0% (10 μL) 500 μL: ± 1.0% (5.0 μL) 100 μL: ± 2.5% (2.5 μL)
	Random Error (per channel)	20 $\mu$ L: $\leq$ 0.80% (0.16 $\mu$ L) 10 $\mu$ L: $\leq$ 1.0% (0.10 $\mu$ L) 2 $\mu$ L: $\leq$ 5.0% (0.10 $\mu$ L) 1 $\mu$ L: $\leq$ 10% (0.10 $\mu$ L)	200 $\mu$ L: $\leq$ 0.40% (0.80 $\mu$ L) 100 $\mu$ L: $\leq$ 0.80% (0.80 $\mu$ L) 20 $\mu$ L: $\leq$ 1.5% (0.30 $\mu$ L) 5 $\mu$ L: $\leq$ 3.5% (0.175 $\mu$ L)	' ' ' '
	Volume Increments	0.02 μL	0.2 μL	1 μL

	VOLUME RANGE	2 – 20 μL	5 – 200 µL
96.	Systematic Error (per channel)	20 μL: ± 1.0% (0.20 μL) 10 μL: ± 1.2% (0.12 μL) 2 μL: ± 6.0% (0.12 μL) 1 μL: ± 12% (0.12 μL)	200 μL: ± 1.0% (2.0 μL) 100 μL: ± 1.0% (1.0 μL) 20 μL: ± 2.0% (0.40 μL) 5 μL: ± 5.0% (0.25 μL)
Liquidator 96	Random Error (per channel)	20 $\mu$ L: $\leq$ 0.80% (0.16 $\mu$ L) 10 $\mu$ L: $\leq$ 1.0% (0.10 $\mu$ L) 2 $\mu$ L: $\leq$ 5.0% (0.10 $\mu$ L)	200 $\mu$ L: $\leq 0.50\%$ (1.00 $\mu$ L) 100 $\mu$ L: $\leq 0.80\%$ (0.80 $\mu$ L) 20 $\mu$ L: $\leq 1.5\%$ (0.30 $\mu$ L)

 $1 \mu L$ :  $\leq 10\% (0.12 \mu L)$ 

0.02 μL

200 μL

0.2 μL

 $5~\mu L: \leq 3.5\% \quad (0.175~\mu L)$ 

	MODEL/MAX. VOLUME	20 μL	300 μL	1250 μL
MicroPro	VOLUME RANGE	2 – 20 μL	5 – 300 µL	50 – 1250 μL
	Systematic Error (per channel)	20 μL: ± 1.0% (0.20 μL) 10 μL: ± 1.5% (0.15 μL) 2 μL: ± 5.0% (0.10 μL)	300 μL: ± 1.0% (3.0 μL) 150 μL: ± 1.5% (2.25 μL) 30 μL: ± 1.5% (0.45 μL) 15 μL: ± 2.0% (0.30 μL)	1250 μL: ± 1.8% (22.5 μL) 625 μL: ± 2.0% (12.5 μL) 125 μL: ± 2.5% (3.13 μL)
	Random Error (per channel)	20 $\mu$ L: $\leq 0.5\%$ (0.10 $\mu$ L) 10 $\mu$ L: $\leq 1.0\%$ (0.10 $\mu$ L) 2 $\mu$ L: $\leq 3.0\%$ (0.060 $\mu$ L)	$300 \ \mu L$ : $\leq 1.0\% \ (3.0 \ \mu L)$ $150 \ \mu L$ : $\leq 1.0\% \ (1.50 \ \mu L)$ $30 \ \mu L$ : $\leq 1.0\% \ (0.30 \ \mu L)$ $15 \ \mu L$ : $\leq 2.0\% \ (0.30 \ \mu L)$	1250 μL: $\leq$ 0.8% (10.0 μL) 625 μL: $\leq$ 1.0% (6.25 μL) 125 μL: $\leq$ 1.5% (1.88 μL)
	Volume Increments	0.01 μL	0.1 μL	1 μL

## Rainin QuickFlow

# Portable Easy-to-Use Aspirator

This compact benchtop aspirator delivers reliable performance while combining exceptional ergonomics with a broad range of features and conveniences.



# Rainin SP+ Pipette Controller

# Easy Control at Your Fingertips

Superior ergonomics, state-of-the-art technology and a unique, extendable head to conquer your larger volume liquid handling needs with speed and ease.

#### Extendable head

Telescoping head for added comfort, particularly when operating in a laminar flow hood.

#### Variable speed, low force

Easily adjustable aspiration/dispense speeds offer greater flexibility with various solutions and media, and protection for sensitive cells.

#### Lightweight, easy to use

The ergonomic and well-balanced design makes the instrument easy to hold and operate.

#### **Powerful battery**

Delivers 15 hours of operation on single charge and up to one hour after just 10 minutes.

#### Wide range of consumables

Rainin SP+ works with all types of 1 mL to 100 mL serological and micro pipettes and can transfer liquid from 0.1 mL to 100 mL.

	CATALOG No.	DESCRIPT	ION			PRICE		
+	30671415	Rainin SF	+ pipette co	ontroller, 0.1–10	0 mL	\$ 693.00		
SP+	17011745	Rainin SF	+ spare filt	er		17.40		
Rainin	30671310	Magnetic	Hang-Up			62.40		
ž	30621247	Pipette ac	Pipette adapter					
	CATALOG No.	SIZE	MAX. VOL.	QUANTITY		PRICE		
es	30749261	1 mL	1.3 mL	200 per pack	Sterilized	\$ 45.80		
pett	30749262	2 mL	2.8 mL	200 per pack	Sterilized	53.80		
<u>-</u>	30749263	5 mL	8 mL	100 per pack	Sterilized	32.20		
gica	30749264	10 mL	13 mL	100 per pack	Sterilized	37.80		
Serological Pipettes	30749265	25 mL	33 mL	50 per pack	Sterilized	37.80		
Se	30749266	50 mL	60 mL	50 per pack	Sterilized	80.80		



Serological pipettes
Color coded with easy to read graduation marks.

mt.com/SR-SPplus

# Rainin AutoRep S

# Manual Repeater Pipette

Flexible volume ranges and a wide selection of syringes ensure precise performance.

The ergonomic AutoRep™ S manual repeater pipette is simple to operate, lightweight and easy on the hands.

Useful when working in the aliquot range from 2 µL to 5 mL, AutoRep S is an economical choice for repeat dispensing up to 49 aliquots without refilling.

The wide selection of Encode™ syringes in sizes from 0.1 mL to 50 mL is available in sterile and non-sterile.

	CATALOG No.	DESCRIPTION			PRICE		
S de	17013008		Manual repeater pipette volume range: 2 µL to 5 mL				
AutoRep S	17013009	Shelf mount	Shelf mount				
Α	17013010	Volume table			40.80		
	CATALOG No.	MAX. VOLUME	QUANTITY				
	17007399	0.1 mL	100 per pkg		\$ 147.20		
	17007400	0.1 mL	100 per pkg	Sterile	207.20		
	17001871	0.5 mL	100 per pkg		147.20		
	17001872	0.5 mL	100 per pkg	Sterile	207.20		
	17013004	1 mL	100 per pkg		147.20		
	17013002	1 mL	100 per pkg	Sterile	210.20		
	17001873	1.25 mL	100 per pkg		147.20		
10	17001874	1.25 mL	100 per pkg	Sterile	207.20		
Encode Syringes	17001877	2.5 mL	100 per pkg		147.20		
Syri	17001878	2.5 mL	100 per pkg	Sterile	207.20		
de	17001883	5 mL	100 per pkg		147.20		
inco	17001884	5 mL	100 per pkg	Sterile	207.20		
_	17013005	10 mL	100 per pkg		156.20		
	17013003	10 mL	100 per pkg	Sterile	210.20		
	17001875	12.5 mL	100 per pkg		147.20		
	17001876	12.5 mL	100 per pkg	Sterile	207.20		
	17001879	25 mL*	50 per pkg		168.60		
	17001880	25 mL*	25 per pkg	Sterile	328.00		
	17001881	50 mL*	25 per pkg		167.20		
	17001882	50 mL*	25 per pkg	Sterile	328.00		





# Rainin Disp-X

# Large Volume Bottle Top Dispensing

Built-in safety features prevent spills and ensure smooth, trouble-free operation.

#### Superior chemical compatibility

High-grade inert materials ensure the purity of your reagents from the first aliquot to the last.

#### Safety first!

Discharge tube and closure cap prevent accidental dispensing or spilling of corrosive liquids.

#### Clear and easy control

Adjustable control knob facilitates volume setting.

#### Flexible options

Telescoping filling tubes accommodate different size bottles.



	CATALOG No.	DESCRIPTION			PRICE
Disp-X	30373527	Disp-X 0.5-5 mL bottle-top dispenser			\$ 522.00
	30373528	Disp-X 1-10 mL bottle-top dispenser			522.00
	30373529	Disp-X 2.5-25 mL bottle-top dispenser			714.00
	30373750	Disp-X 5-50 mL bottle-top dispenser			714.00
Accessories	30373751	Disp-X discharge tube	5 or 10 mL		\$ 86.40
	30373752	Disp-X discharge tube	25 or 50 mL		86.40
	30373753	Disp-X filling tube	5 or 10 mL	125-240 mm length	34.60
	30373754	Disp-X filling tube	5 or 10 mL	250-480 mm length	34.60
	30373755	Disp-X filling tube	25 or 50 mL	170-330 mm length	34.60
	30373756	Disp-X filling tube	25 or 50 mL	250-480 mm length	34.60
	30373757	Disp-X filling valve	5 or 10 mL		86.40
	30373758	Disp-X filling valve	25 or 50 mL		85.80
	30373759	Disp-X discharge valve	5 or 10 mL		86.40
	30373760	Disp-X discharge valve	25 or 50 mL		86.80

# **Rainin**Service

# Our Lab or Yours - You Choose.

#### Mail-in

- Service when you need it no waiting to schedule!
- Fast turnground
- ISO 17025 accredited by A2LA
- · Tightly-controlled environment
- 6- and 7-place METTLER TOLEDO balances
- METTLER TOLEDO MCP multichannel balances
   all channels tested
- · Free inbound shipping
- Free return kits (4+ pipettes)

# **Onsite**

- · We come to you!
- · Perfect for labs with chain-of-custody concerns
- ISO 17025 accredited by A2LA
- METTLER TOLEDO balances
- Minimal downtime
- Pipettes are calibrated in the environment where they are used
- Experienced Rainin onsite technicians available nationwide

Mail-in calibration services now available on the ShopRainin e-store. Everything you need to manage mail-in pipette calibrations now in one convenient location.

▶ ShopRainin.com/service



# Why choose Rainin Service? Flexibility, convenience and highly trained technicians!

#### Choice

Our service offering includes onsite or mail-in calibration service and a wide range of plans to choose from.

#### Trust

Our highly trained service experts have years of experience repairing and calibrating many different brands of pipettes.

We only use Rainin parts on Rainin pipettes.

#### Confidence

Rainin pipette service has the knowledge, tools and certifications to deliver expert performance every time.

Every calibration includes meticulous cleaning of each component, part replacement as needed, and preventive maintenance.

Can't decide? Try our Service Finder to see which plans fit your lab's needs.

▶ ShopRainin.com/service-finder



# Mail-in Pipette Calibration on ShopRainin.com

# Your One Stop Shop for Mail-in Calibration Services

As the world's largest pipette service provider, Rainin makes it easy to get fast reliable service and calibration, guaranteed!

- Quick turnaround, with four ISO 17025accredited service labs across North America
- We service all brands of pipettes most with manufacturers' original parts (Rainin®, Eppendorf®, Gilson®, Integra®, Thermo Fisher®, Sartorius®)
- We service all pipettes single channels, multichannels and all Rainin benchtop pipettors
- · Free inbound shipping
  - 4 or more, request a free return kit
  - 3 or fewer, or in a rush, pack them yourself and print a complimentary inbound shipping label

Order Mail-in calibration service online:

▶ ShopRainin.com/service

#### Service Plans

#### All plans include: Standard Plan **Advanced Plan Accredited Plan** Accredited+ Plan · Full preventive maintenance (PM) with seal replacement or restoration Labs with basic When prior For complete ISO 17025 ISO 8655 and 17025 · Liquid leak test on all service needs performance matters documentation compliance nozzles · Adjustment to manufacturer's Pharma tolerances GLP/GMP, R&D · Calibration certificate with Clinical Regulated Labs Regulated Labs detailed results esearch Labs Research Labs · Calibration label 4 Test weighings @ 4 Test weighings @ 10 Test weighings @ As Found N/A 10, 100% 10, 50, 100% 10, 50, 100% 4 Test weighings @ 4 Test weighings @ 10 Test weighings @ 10 Test weighings @ As Returned 10, 100% 10, 100% 10, 50, 100% 10, 50, 100% Preventive Maintenance Yes Yes Yes Yes **Calibration Certificate** Non-accredited Non-accredited Accredited Accredited Compliance ISO 8655-7 Annex A ISO 8655-7 Annex A ISO 8655-7:2022 ISO 8655-6:2022

# Mail-in Ordering Information

# As Found N/A As Returned 4 Test weighings @ 10, 100%

- Full preventive maintenance (PM) with leak test and seal replacement or restoration
- · Adjustment to manufacturer's tolerances
- Calibration certificate with detailed results and historical trending

CATALOG No.	DESCRIPTION	PRICE
17800031	Single channel	\$ 53.40
17800096	Multichannel All — 8 channels	143.00
17800131	Multichannel All — 12 channels	171.20
17800153	Liquidator 96	2,475.00
30451098	BenchSmart 96 Head	2,475.00
30848615	MicroPro	2,475.00

Accredited Plan	
As Found	4 Test weighings @ 10, 50, 100%
As Returned	10 Test weighings @ 10, 50, 100%

This plan includes all services provided under the Advanced Plan. **PLUS:** 

- As found A2LA calibration results and failure confirmation
- A2LA calibration certificate incl. uncertainty of measurement (ISO 17025 compliant), historical trending and calibration label

CATALOG No.	DESCRIPTION	PRICE
17800048	Single channel	\$ 124.60
17800113	Multichannel All — 8 channels	378.00
17800148	Multichannel All — 12 channels	423.00
30415512	Liquidator 96	3,885.00
30451100	BenchSmart 96 Head	3,885.00
30848617	MicroPro	3,885.00

# Advanced Plan As Found 4 Test weighings @ 10, 100% As Returned 4 Test weighings @ 10, 100%

This plan includes all services provided under the Standard Plan.

#### PLUS:

 As found calibration results, including failure confirmation

CATALOG No.	DESCRIPTION	PRICE
17800040	Single channel	\$ 64.00
17800105	Multichannel All — 8 channels	197.20
17800140	Multichannel All — 12 channels	241.00
17800154	Liquidator 96	2,920.00
30451099	BenchSmart 96 Head	2,920.00
30848616	MicroPro	2,920.00

Accredited+ Plan	
As Found	10 Test weighings @ 10, 50, 100%
As Returned	10 Test weighings @ 10, 50, 100%

This plan includes all services provided under the Accredited Plan.

#### PLUS:

 Full compliance to ISO 8655:2022 part 6 and meets all quality standards

CATALOG No.	DESCRIPTION	PRICE
30829339	Single channel	\$ 179.20
30829341	Multichannel All — 8 channels	571.00
30829342	Multichannel All — 12 channels	679.00

Optional Services		
DNA/RNase dec	contamination	
CATALOG No.	DESCRIPTION	PRICE
17800171	Single channel	\$ 30.80
17800401	Multichannel	61.40
30932020	96-Channel	184.00

# **Onsite Pipette Calibration**


## **Expert Calibration in Your Lab**

Rainin makes it easy to receive state-of-the-art ISO 17025-accredited pipette service right in your lab. We service all brands.

Regardless of plan, all pipettes we service receive:

- Full preventive maintenance (PM) with seal replacement or restoration
- Liquid leak test on all nozzles
- Calibrated to ISO 8655:2022 tolerances
- Calibration certificates with detailed results (except Basic Plan)
- Calibration label

#### Service Plans

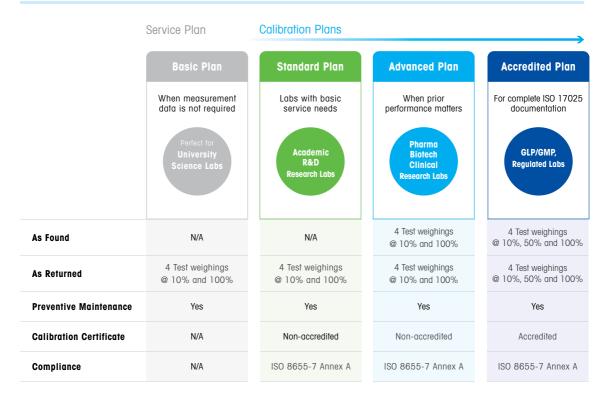

To schedule onsite service, call 800.472.4646 (choose option "2" or Sales).

# Onsite Ordering Information

Basic Plan	
As Found	N/A
As Returned	4 Test weighings @ 10, 100%

• Full preventive maintenance (PM) with leak test and seal replacement or restoration

CATALOG No.	DESCRIPTION	PRICE*
30445525	Single channel	\$ 39.80
30445526	Multichannel All — 8 channels	75.40
30445527	Multichannel All – 12 channels	93.60

	Standard Plan
As Found	N/A
As Returned	4 Test weighings @ 10, 100%
Full preventive maintenance (PM) with leak test and seal replacement or restoration	

· Calibration certificate with detailed

measurement data

CATALOG No.	DESCRIPTION	PRICE*
17800347	Single channel	\$ 53.40
17800292	Multichannel All — 8 channels	146.60
17800256	Multichannel All — 12 channels	175.20

	Advanced Plan
As Found	4 Test weighings @ 10, 100%
As Returned	4 Test weighings @ 10, 100%

This plan includes all services provided under the Standard Plan.

#### PHIS

As found calibration results – including failure confirmation

CATALOG No.	DESCRIPTION	PRICE*
17800320	Single channel	\$ 65.20
17800272	Multichannel All — 8 channels	201.20
17800236	Multichannel All – 12 channels	246.00

	Accredited Plan
As Found	4 Test weighings @ 10, 50, 100%
As Returned	4 Test weighings @ 10, 50, 100%

This plan includes all services provided under the Advanced Plan.

#### PLUS

ISO 17025 calibration certificate, including uncertainty of measurement (certified by A2LA), historical trending and calibration label

CATALOG No.	DESCRIPTION	PRICE*
30451900	Single channel	\$ 114.80
30451876	Multichannel All — 8 channels	317.00
30451852	Multichannel All – 12 channels	399.00

 $<sup>^{</sup>st}$  Travel fees may apply. Contact your sales representative.

# **Expert Calibration and Repair**

# Buy Vouchers Now, Use Them When You Need Them

Purchase service vouchers when you buy your pipettes or when you have funding. Use them later — whenever you're ready for calibration!

- · Valid for two years
- Use with any brand of pipette
- Eliminate separate purchase orders
- Purchase with your pipettes and tips
- Valid for mail-in service only

	CATALOG No.	DESCRIPTION	PRICE						
	17014448	Standard Single channel service voucher	\$ 53.40						
	17014453	Standard Multichannel All 8-channel service voucher	143.00						
	17014456	Standard Multichannel All 12-channel service voucher	171.20						
	30415513	Standard Liquidator 96 – 96 channels, includes box & round-trip UPS Ground shipping	2,475.00						
	30456199	Standard BenchSmart 96 Head — 96 channels, includes round-trip UPS Ground shipping							
က်	17014449	Advanced Single channel service voucher							
Vouchers	17014454	Advanced Multichannel All 8-channel service voucher	197.20						
	17014457	Advanced Multichannel All 12-channel service voucher	241.00						
Service	30415514	Advanced Liquidator 96 – 96 channels, includes box & round-trip UPS Ground shipping	2,920.00						
Se	30456200	Advanced BenchSmart 96 Head — 96 channels, includes round-trip UPS Ground shipping	2,920.00						
	17014450	Accredited Single channel service voucher	124.60						
	17014455	Accredited Multichannel All 8-channel service voucher	378.00						
	17014458	Accredited Multichannel All 12-channel service voucher	423.00						
	30415517	Accredited Liquidator 96 – 96 channels, includes box & round-trip UPS Ground shipping	3,885.00						
	30456201	Accredited BenchSmart 96 Head — 96 channels, includes round-trip UPS Ground shipping	3,885.00						



## At Your Service!


# We're Here to Help



#### Pipette Inventory Assessment

Time for a pipette review? Contact your local Rainin representative to request an inventory assessment. You'll receive a complimentary summary report that includes quantity, age and ergonomic evaluation. Optimize your pipette inventory to your current laboratory workflow, eliminate redundancy and minimize your risk of using out-of-calibration pipettes.

www.mt.com/pipette-inventory

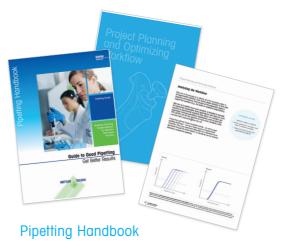

This guide provides an in-depth overview of pipettes, tips, pipetting techniques and pipetting accuracy for both new and experienced users.

Pipetting is probably the most frequently practiced activity in life science labs today. Understanding the basics of good pipetting practice is important to the success and reproducibility of any experiment. No matter how precise a pipette is, the skill and knowledge of the user ultimately determines the accuracy and reliability of their results.

www.mt.com/pipetting-handbook



- Contact Your Sales Representative 800 4 RAININ (800-472-4646) or is@rainin.com
- Contact Technical Support or Service
   Reach our technical support team at 800 4 Rainin (800-472-4646, choose option "3") or by email at tech.support@rainin.com

# **Good Pipetting Practice**

# Better Results Through Better Understanding

The performance of any instrument is improved in the hands of a skilled operator. Through GPP<sup>™</sup>, labs can reduce errors and improve the reliability and repeatability of their results. The key: better understanding the factors that affect experiments and learning application-specific techniques.

You can improve data quality with Good Pipetting Practice™ – METTLER TOLEDO's comprehensive, systematic approach to maximizing pipetting accuracy and repeatability.



injury through better ergonomics.

## **GPP Toolbox**

# A Wide Range of Tools



#### **GPP Seminars**

We offer liquid handling-related seminars to meet the needs of every stakeholder in your organization, from bench scientists and lab managers to quality control, health and safety. We hold short and long format trainings online or in in your lab and can even measure each user's pipetting proficiency before and after training.

www.mt.com/gpp



## **Technique Posters**

From good techniques to pipette cleaning to good ergonomics, our free technique posters reinforce pipetting practices that will help ensure more consistent and accurate results.

www.mt.com/rainin-posters





Scan the QR code to discover what GPP can do for your lab!

## **Reduce Plastic Waste with Rainin**

# **Greener Tip Rack Solutions**

Rainin offers a number of waste-saving pipette tip rack options to help meet your sustainability objectives. We've taken a three-pronged approach to reducing plastic rack waste:

- Reducing the amount of plastic particularly #5 polypropylene that, when recycled, downcycles into lower-quality products (e.g., doormats, crates, etc.).
- Replacing polypropylene with PET, which is more widely recycled. Recycled #1 PET can be
  made into the same or similar products (water bottles, TerraRacks, etc.).
- · Avoiding plastic components, instead using natural, compostable materials like plant fiber.

Plastic savings shown below are calculated using our current standard 1000  $\mu$ L polypropylene tip rack as the baseline. Components used for the calculation include the rack base, lid, tip deck, sleeve (refills) and any plastic wrapping – tips are not included. Actual plastic savings percentages vary slightly by rack size (10  $\mu$ L racks are smaller than 1000  $\mu$ L racks). Our current standard rack, which we introduced in 2018, uses 20% less polypropylene than the rack it replaced, so the amount of plastic saved shown below is on top of the already-reduced plastic rack.



#### **Bulk Tips**

The most sustainable rack is the rack you don't use, which makes manually refilling spent racks and wafers with BioClean Ultra tips in bulk a reasonable option for some labs. See page 55.




#### EarthRack

EarthRack eliminates plastic waste almost entirely. Other than the thin polypropylene tip deck, EarthRack is made entirely from compostable, rapidly renewable plant fiber. EarthRack uses up to 82% less plastic than the Rainin standard rack. See page 48.



#### Green-Pak

Green-Pak and Green-Pak SpaceSaver refill systems for Rainin standard racks reduce plastic waste by up to 80%. Green-Pak's SpaceSaver sleeve and Green-Pak sealed refill packs are both made from highly recyclable PET. See page 52.

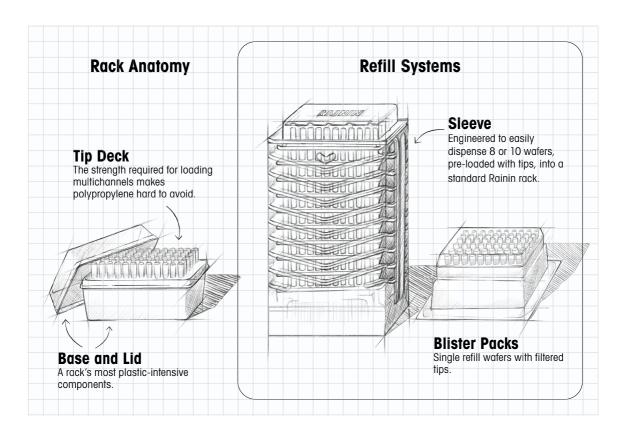



#### StableStak

Rainin StableStak multilayer system is designed for high-production multichannel work. With five layers mounted on a sturdy base and 192 tips per layer, each StableStak is equivalent to 10 96-tip racks. See page 54.



#### TerraRack

Not only composed of 50% less plastic than the Rainin standard rack, TerraRack is made almost entirely from PET – the same plastic as in water bottles and the most commonly recycled plastic worldwide. Easy to disassemble and compress, TerraRacks take up considerably less volume as waste. See page 50.



#### Rainin Standard Rack

Like most pipette tip racks, Rainin standard racks are 100% polypropylene (#5), which is prized for strength and durability. This rack was reengineered in 2018 to reduce the amount of polypropylene it contains by 20%. See page 56.

### Rainin EarthRack

# Compostable, Biodegradable Tip Rack

EarthRack is made primarily from refined cellulose that is easy to compost and is biodegradable. Because the small, recyclable polypropylene tip deck is the only plastic component, EarthRack significantly reduces plastic use and waste.

Available in 20  $\mu$ L, 200  $\mu$ L and 1000  $\mu$ L. Rainin LTS LiteTouch System only.





Scan the QR code for a timelapse video of EarthRack biodegrading.





#### Earth friendly

EarthRack's lid and base are composed of a biodegradable bamboo and sugar cane cellulose material that can be composted at municipal facilities. If deposited in a landfill, they will biodegrade.



### Reduces plastic waste

With up to 75% less plastic than conventional racks, adopting EarthRack can help organizations meet their sustainability objectives.



#### Less carbon

When composted, incinerated or deposited into a landfill, EarthRack's cellulose components contribute no net increase in atmospheric CO<sub>2</sub>. Carbon absorbed during initial plant growth is simply returned to the atmosphere.



#### **Specifications**

BioClean	Green	Test	Limits	Ωf	Detection	(LODs)

Human DNA	<0.32 pg
Bacterial DNA	<0.032 pg
Plant (chloroplast) DNA	<0.032 pg
DNase (three tests)	
Gel electrophoresis	≤10 <sup>-7</sup> KU/µL
• Kunitz Assay	≤10 <sup>-4</sup> KU/µL
• qPCR	≤10 <sup>-5</sup> KU/µL
RNase (three tests)	
Gel electrophoresis	≤10 <sup>-5</sup> KU/µL
• Kunitz Assay	≤10 <sup>-5</sup> KU/µL
• RT-qPCR	≤10 <sup>-5</sup> KU/µL
Endotoxin (LPS)	≤0.01 EU/mL
Protein	<2 ng/tip
Protease	<500 ng/mL
PCR inhibitors	None detected

# EarthRack is **Bio**Clean Green™

Rainin EarthRack products are certified to be BioClean Green - a comprehensive testing process and specifications developed specifically for cellulose-based tip racks.

Because EarthRack is derived from natural, plant-based materials, BioClean Green employs a broader, more intrusive quality testing procedure to ensure the absence of biological contaminants. For example, the use of cellulose requires an additional test for chloroplast DNA. Testing for DNase and RNase has been expanded to include three separate tests for each nuclease for complete assurance of tip performance.

				\	'OL	QTY		FEATUR	RES
		CATALOG No.	DESCRIPTION	M44 VO/	PACKS.	TPSCART	STEP.	GENTE	PRICE
LTS	Š	30807966	ER LTS 20µL F 960G/10	20 μ	L 10	960	✓	✓	\$ 180.40
	EarthRack	30807967	ER LTS 200µL F 960G/10	200 μ	L 10	960	✓	✓	180.40
	Ear	30808038	ER LTS 1000µL F 768G/8	1000 μ	L 8	768	✓	✓	180.40

mt.com/SR-BCG



## **TerraRack**

# The Most Eco-Friendly Rack Around

Less weight, less waste, 100% recyclable

TerraRack™ is a revolutionary concept in tip racks. Sturdy as conventional racks, yet made with less than half the plastic of standard racks and completely recyclable. In fact, the hinged TerraRack shell is made from PET, which is easily recycled.

For convenience and complete cleanliness, TerraRacks with BioClean Ultra tips come sterilized (they are not autoclavable).

Many volume sizes available in both LTS [115] and universal unv formats.







#### Less waste

With 50% less plastic than conventional racks, TerraRacks are also extremely compressible, so racks occupy significantly less space in a waste or storage container.



#### 100% recyclable

In addition to being completely recyclable, much of each TerraRack is composed of PET, which is easily recycled and aggressively recovered from the plastic waste stream.



#### Non slip accessory

Each TerraBase provides secure footing on slippery bench surfaces for an enhanced pipetting experience. Durable construction provides high-impact resistance.

mt.com/SR-BCU





# **Eco-Friendly Refills**

# Less Waste, More Space

Space-saving rack refill options enhance your workflow and minimize plastic waste. Order empty racks if you don't already have them. Many volume sizes available in both LTS us and universal workflow formats.





Green-Pak®						QT	Υ	FE	ATURE	S
Individually seale	ed re	fills — idea	Il for filter tip users	Max VOLII.	JME	TIPS CARTON	STEE	2	24	
				\\ \\ \\ \\ \\ \\ \\ \\ \\ \\ \\ \\ \\	/ c		<del>7</del> / 4	SHAE	FILTE	& /
		CATALOG No.	DESCRIPTION	71	\$		/ 5	. / 55	/ 2	PRICE
LTS		30389270	GP LTS 20µL 960A/10	20 μL¹	10	960				\$ 81.60
		30389275	GP LTS 20µL S 960A/10	20 μL¹	10	960	<b>✓</b>			105.80
		30389274	GP LTS 20µL F 960A/10	20 μL¹	10	960	<b>✓</b>		<b>√</b>	162.80
	Pak	30389277	GP LTS 250µL 960A/10	200 μL <sup>2</sup>	10	960				81.60
	Green-Pak	30389278	GP LTS 250µL S 960A/10	200 µL²	10	960	<b>✓</b>			105.80
	Gre	30389276	GP LTS 200µL F 960A/10	200 μL <sup>2</sup>	10	960	<b>✓</b>		<b>√</b>	162.80
		30389271	GP LTS 1000µL 768A/8	1000 μL	8	768				81.60
		30389273	GP LTS 1000µL S 768A/8	1000 μL	8	768	<b>✓</b>			105.80
		30389272	GP LTS 1000µL F 768A/8	1000 µL	8	768	✓		$\checkmark$	162.80
UNV		30389258	GP UNV 10µL 960A/10	10 µL	10	960				\$ 81.60
		30389265	GP UNV 10µL S 960A/10	10 µL	10	960	✓			105.80
		30389263	GP UNV 10µL F 960A/10	10 µL	10	960	<b>✓</b>		✓	162.80
		30389264	GP UNV 10µL FG 960A/10	10 µL	10	960	✓	<b>✓</b>	✓	162.80
	¥	30389267	GP UNV 20µL F 960A/10	20 µL	10	960	✓		✓	162.80
	Green-Pak	30389262	GP UNV 100µL F 960A/10	100 μL	10	960	✓		✓	162.80
3038926		30389266	GP UNV 200µL F 960A/10	200 μL	10	960	✓		✓	162.80
	9	30389268	GP UNV 250µL 960A/10	250 µL	10	960				81.60
		30389269	GP UNV 250µL S 960A/10	250 µL	10	960	✓			105.80
		30389259	GP UNV 1000µL 768A/8	1000 µL	8	768				81.60
		30389261	GP UNV 1000μL S 768A/8	1000 µL	8	768	✓			105.80
		30389260	GP UNV 1000μL F 768A/8	1000 µL	8	768	✓		✓	162.80





#### Fast

Refill a rack in seconds

#### Clean

Protects tips from contamination

#### Convenient

Dispensing sleeve is easy to use

Green-Pak Spac		• •	)	OL	QT		FE	ATURE	S
Each sleeve conta	ains 8 or 19	O refills  DESCRIPTION	MAY VOLIM.	RACE	TIPS/CARTON	STED	SHAE	FILTES	PRICE
LTS	30389291	GPS LTS 20µL 960A/10	20 μL <sup>1</sup>	10	960				\$ 81.60
LIS	30389297	GPS LTS 20µL S 960A/10	20 μL <sup>1</sup>	10	960	/			105.80
-	30389299	GPS LTS 250µL 960A/10	200 μL <sup>2</sup>	10	960				81.60
Sav	30389301	GPS LTS 250µL S 960A/10	200 μL <sup>2</sup>	10	960	<b>✓</b>			105.80
SpaceSaver	30389303	GPS LTS 300µL 768A/8	300 µL	8	768				81.60
Spi	30389304	GPS LTS 300µL S 768A/8	300 µL	8	768	<b>✓</b>			105.80
	30389292	GPS LTS 1000µL 768A/8	1000 μL	8	768				81.60
	30389294	GPS-LTS 1000µL S 768A/8	1000 μL	8	768	<b>✓</b>			105.80
iks	30389354	Tip Racks GPR LTS 20µL 0/10	20 µL	10					\$ 41.60
Rac	30389353	Tip Racks GPR XXX 250µL 0/10	200 μL	10					41.60
Empty Racks	30397667	Tip Racks GPR XXX 300µL 0/8	300 µL	8					41.60
置	30389351	Tip Racks GPR XXX 1000µL 0/8	1000 µL	8					41.60
UNV	30389283	GPS UNV 10µL G 960A/10	10 µL	10	960		<b>√</b>		\$ 81.60
	30389285	GPS UNV 10µL GS 960A/10	10 μL	10	960	<b>✓</b>	✓		105.80
le l	30389287	GPS UNV 250µL 960A/10	250 μL	10	960				81.60
Sav	30389289	GPS UNV 250µL S 960A/10	250 μL	10	960	<b>✓</b>			105.80
SpaceSaver	30397665	GPS UNV 300µL 768A/8	300 µL	8	768				81.60
ďs	30397666	GPS UNV 300µL S 768A/8	300 µL	8	768	✓			105.80
	30389279	GPS UNV 1000µL 768A/8	1000 μL	8	768				81.60
	30389281	GPS UNV 1000µL S 768A/8	1000 μL	8	768	<b>√</b>			105.80
cks	30389350	Tip Racks GPR UNV 10µL 0/10	10 μL	10					\$ 41.60
Ra	30389353	Tip Racks GPR XXX 250µL 0/10	250 μL	10					41.60
Empty Racks	30397667	Tip Racks GPR XXX 300µL 0/8	300 µL	8					41.60
盂	30389351	Tip Racks GPR XXX 1000µL 0/8	1000 μL	8					41.60



# **Double-96 Racked Tips**

## StableStak and StableRak

Ideal for high production environments. Perfect for multichannel pipetting! Available in both LTS us and universal w formats.

StableStak and StableRak are available in standard, sterilized and filter (StableRak only) formats. They keep the tips safe from debris, dust and contamination. Their sturdy construction using non-slip racks allows fast and easy loading.

StableRak™			StableStak™							
Individual double	-96	racks	Stacked double-96 racks							
with 192 tips per	lave	er.	with 192 tips per layer.	VO		QT	Υ	FEA	TURES	
		CATALOG No.	DESCRIPTION	Max ! OLUM.	RACK	7 IPS/CARTON	STEPH	FILTER	SHAFTGARD	PRICE
LTS		17005862	SR LTS 20µL 960A/5	20 μL¹	5	960				\$ 81.40
		17005861	SR LTS 20µL S 960A/5	20 μL¹	5	960	✓			105.80
		17005860	SR LTS 20µL F 960A/5	20 μL¹	5	960	<b>✓</b>	<b>√</b>		162.80
		17005864	SR LTS 250µL 960A/5	200 μL²	5	960				81.40
Section 1		17005863	SR LTS 250µL S 960A/5	200 μL <sup>2</sup>	5	960	✓			105.80
		17005859	SR LTS 200µL F 960A/5	200 μL²	5	960	✓	✓		162.80
	¥	17005867	SR LTS 300µL 768A/4	300 μL	4	768				81.40
	StableRak	17005866	SR LTS 300µL S 768A/4	300 µL	4	768	✓			105.80
	Sta	17005865	SR LTS 300µL F 768A/4	300 μL	4	768	<b>✓</b>	<b>√</b>		162.80
		17007083	SR LTS 1000µL 768A/4	1000 μL	4	768				81.40
		17007082	SR LTS 1000µL S 768A/4	1000 μL	4	768	<b>✓</b>			105.80
		17007081	SR LTS 1000µL F 768A/4	1000 μL	4	768	✓	✓		162.80
		17007086	SR LTS 1200µL 768A/4	1200 μL	4	768				81.40
		17007085	SR LTS 1200µL S 768A/4	1200 μL	4	768	✓			105.80
		17007084	SR LTS 1200µL F 768A/4	1200 µL	4	768	✓	✓		162.80
		17005873	SS LTS 20µL 960A/5	20 μL¹	5	960				\$ 81.40
		17005872	SS LTS 20µL S 960A/5	20 μL¹	5	960	✓			105.80
	품	17005875	SS LTS 250µL 960A/5	200 μL <sup>2</sup>	5	960				81.40
	StableStak	17005874	SS LTS 250µL S 960A/5	200 μL²	5	960	✓			105.80
	g	17005877	SS LTS 300µL 768A/4	300 µL	4	768				81.40
	S	17005876	SS LTS 300µL S 768A/4	300 µL	4	768	<b>✓</b>			105.80
		17007090	SS LTS 1000µL 768A/4	1000 μL	4	768				81.40
		17007089	SS LTS 1000µL S 768A/4	1000 μL	4	768	✓			105.80
UNV		17005857	SR UNV 300µL 768/4	300 µL	4	768				\$ 81.40

# **Bulk Tips**

# The Economic Solution

Ideal for academia and other research laboratories. Available in both LTS us and universal w formats.

				Max V0/	PAC. JUNE	The Share of the state of the s	7	FILE	EATURE ON SHAPE	WIDE ONE OF	PRICE
LTS		17001128	RC LTS 20µL 1000A/1	20 μL¹		1000					\$ 62.60
		17001118	RC LTS 250µL 1000A/1	200 μL <sup>2</sup>		1000					62.60
		17001132	RC LTS 300µL 1000A/1	300 µL		1000					62.60
		17001129	RC LTS 1000µL 1000A/1	1000 μL		1000					62.60
	ips	17006324	RC LTS 1200µL 1000A/1	1200 µL		1000					62.60
	Bulk Tips	17001130	RC LTS 2000µL 1000A/1	2000 μL		1000					127.00
	B	17001133	RC LTS 5000µL 1000A/1	5000 μL		1000					142.00
		17001119	RC LTS 10mL 200A/1	10 mL		200					82.40
		17005940*	RC LTS 10mL S 75A/*	10 mL		75	<b>√</b>				133.40
		17001131	RC LTS 20mL 100A/1	20 mL		100					82.40
		17005941*	RC LTS 20mL S 50A/*	20 mL		50	$\checkmark$				98.40
UNV		17001124	RC UNV 10µL 1000A/10	10 µL		1000					\$ 62.60
		17001123	RC UNV 10µL G 1000A/10	10 µL		1000			✓		62.60
		17001116	RC UNV 250μL 1000A/1	250 μL		1000					62.60
	ဟ	17001127	RC UNV 250µL W 1000A/1	250 μL		1000				✓	62.60
	Tips	17001121	RC UNV 1000µL 1000A/1	1000 μL		1000					62.60
	Bulk	17001125	RC UNV 2000µL 1000A/1	2000 μL		1000					127.00
		17001126	RC UNV 2500µL 1000A/1	2500 μL		1000					127.00
		17001117	RC UNV 5000µL 1000A/1	5000 μL		1000					143.20
		17001122	RC UNV 10mL 200A/1	10 mL		200					82.40
		17005939*	RC UNV 10mL S 75A/*	10 mL		75	✓				135.00
		*individually v	vrapped Wide-O = Wide Orifice	<sup>1</sup> Fits 2, 10 (	and 20	μL pipe	ettes	² Fits 5	50, 10	0 and 2	200 µL pipettes



RAININ

DAININ

LTS Sterilia

PRICE

\$81.40

105.80

# **Racked Tips**

# SBS Standard Tip Packaging

#### Removable-cover Racks

Removable-cover racks – easy, practical, convenient packaging. Many volume sizes available in both LTS [15] and universal [10] formats.

- 20% less plastic than the original RT rack!
- Locking latch
- Removable hinged lid
- BioClean Ultra-certified
- Use with Green-Pak and SpaceSaver refills

•	Available with	Low	Retention	tips	(see page	58)

30389239

30389253

30389255

30389254

30389211

30389217

30389215


30389220

30389212

30389230

30389234

30389257

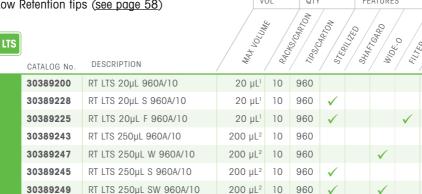



		'
	30389231	RT LTS 1200µL F 768A/8
	30389236	RT LTS 2000µL 480A/8
<sup>1</sup> Fits 2, 10 and	30389238	RT LTS 2000µL S 480A/8
20 µL pipettes	30389237	RT LTS 2000µL F 480A/8
Fits 50, 100 and	30389256	RT LTS 5000µL 192A/8
200 µL pipettes	20200257	DT LTC 5000ul C 1024/9

QTY

VOL

FEATURES



#### Fast loading

Robust secure racks make tip loading easy

#### **Autoclayable**

15 psi, 121°C, 20 minutes


#### Format options

Standard, sterilized and filter versions available

#### **Protective**

Safe from dirt and contamination

mt.com/SR-BCU

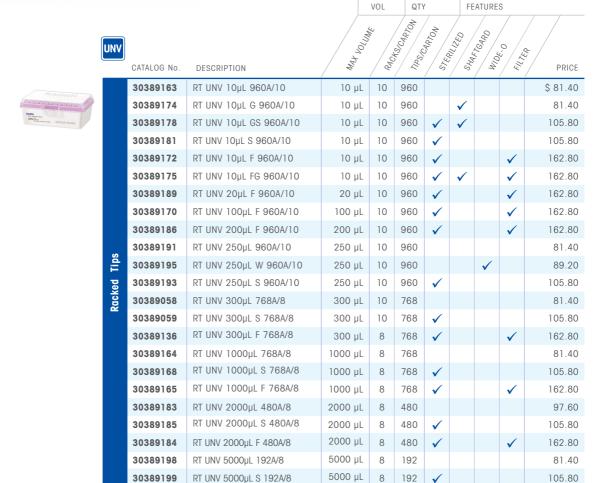

Wide-O = Wide Orifice



# **Low Retention Tips**

# Stop Wasting Valuable Sample!

Available in both LTS and universal works

Superhydrophobic Rainin LR low retention tips minimize residual volume, particularly when transferring viscous samples and liquids with low surface tension. A great option for pipetting valuable molecular biology reagents and expensive samples.

Available in bott	1 LI	s 📧 ana t	iniversal www formats.	V	OL	QT'	1	FE	ATURE	S	
	LTS			MAY VOLLIN.	3M.	TIPS/CARTON	WOJARJON FE	WIG	0.50	FILTE	<u> </u>
		CATALOG No.	DESCRIPTION	\ \ \ \ \ \ \ \ \ \ \ \ \ \ \ \ \ \ \	\ \ \ \ \ \ \ \ \ \ \ \ \ \ \ \ \ \ \		/ 5	7	/ 3	/ 12	PRICE
earner /		30389227	RT LTS 20µL L 960A/10	20 μL <sup>1</sup>	10	960			✓		\$ 105.80
Management of Ma		30389229	RT LTS 20µL LS 960A/10	20 μL <sup>1</sup>	10	960	✓		✓		115.40
		30389226	RT LTS 20µL FL 960A/10	20 μL <sup>1</sup>	10	960	$\checkmark$		✓	<b>√</b>	177.00
		30389244	RT LTS 250µL L 960A/10	200 μL <sup>2</sup>	10	960			✓		105.80
		30389248	RT LTS 250µL LW 960A/10	200 μL <sup>2</sup>	10	960		✓	✓		105.80
		30389246	RT LTS 250µL LS 960A/10	200 μL <sup>2</sup>	10	960	✓		✓		115.40
		30389250	RT LTS 250µL LSW 960A/10	200 μL <sup>2</sup>	10	960	$\checkmark$	✓	✓		142.00
		30389240	RT LTS 200µL FL 960A/10	200 μL <sup>2</sup>	10	960	✓		✓	✓	177.00
¹ Fits 2, 10 and	10	30389241	RT LTS 200µL FLW 960A/10	200 μL <sup>2</sup>	10	960	$\checkmark$	✓	✓	<b>√</b>	177.00
20 μL pipettes	Ξp	30389160	RT LTS 300µL L 768A/8	300 µL	8	768			✓		105.80
² Fits 50, 100 and	Racked Tips	30389161	RT LTS 300µL LS 768A/8	300 µL	8	768	$\checkmark$		✓		115.40
200 μL pipettes	Rac	30389162	RT LTS 300µL FL 768A/8	300 µL	8	768	✓		✓	✓	177.00
		30389214	RT LTS 1000µL L 768A/8	1000 μL	8	768			✓		105.80
		30389216	RT LTS 1000µL LS 768A/8	1000 μL	8	768	$\checkmark$		✓		115.40
		30389213	RT LTS 1000µL FL 768A/8	1000 μL	8	768	$\checkmark$		✓	<b>√</b>	177.00
		30389219	RT LTS 1000µL LW 768A/8	1000 μL	8	768		✓	✓		105.80
		30389221	RT LTS 1000µL LSW 768A/8	1000 μL	8	768	$\checkmark$	✓	✓		142.00
		30389218	RT LTS 1000µL FLW 768A/8	1000 μL	8	768	$\checkmark$	✓	✓	✓	177.00
		30389233	RT LTS 1200µL L 768A/8	1200 μL	8	768			✓		105.80
		30389235	RT LTS 1200µL LS 768A/8	1200 µL	8	768	$\checkmark$		✓		115.40
and the state of t		30389232	RT LTS 1200µL FL 768A/8	1200 µL	8	768	$\checkmark$		✓	✓	177.00
		30389296	GPS LTS 20µL L 960A/10	20 μL¹	10	960			✓		\$ 105.80
		30389298	GPS LTS 20µL LS 960A/10	20 μL <sup>1</sup>	10	960	$\checkmark$		✓		115.20
	SpaceSaver	30389300	GPS LTS 250µL L 960A/10	200 μL <sup>2</sup>	10	960			<b>√</b>		105.80
Doob book and in children	ceS	30389302	GPS LTS 250µL LS 960A/10	200 μL <sup>2</sup>	10	960	$\checkmark$		✓		115.20
Rack base not included with SpaceSaver.	Spa	30398985	GPS LTS 300µL LS 768A/8	300 µL	8	768	$\checkmark$		✓		105.80
See page 53 for		30389293	GPS LTS 1000µL L 768A/8	1000 μL	8	768			✓		105.80
empty racks.		30389295	GPS LTS 1000µL LS 768A/8	1000 μL	8	768	$\checkmark$		✓		115.20
		17014340	RC LTS 20µL L 1000A/1	20 μL		1000			<b>√</b>		\$ 78.60
	Bulk	17014341	RC LTS 250µL L 1000A/1	200 μL		1000			✓		78.60
		17014342	RC LTS 1000µL L 1000A/1	1000 μL		1000			✓		78.60

#### mt.com/SR-BCU

				V		QTY		FE/	ATURES	3		
					4	100	~		/	/	/	///
	UNV			/ 10/	· /,		24PZ		CARD.	0/	REY	÷ /
		CATALOG No.	DESCRIPTION	MAY VOLUM	PA,	7705/CAR7ON	STED	SHAFE	W.D.	0.30 <b>MO</b>	FILTES	PRICE
RAININ		30389180	RT UNV 10µL L 960A/10	10 μL	10	960				✓		\$ 105.80
NAME OF THE PARTY		30389182	RT UNV 10µL LS 960A/10	10 µL	10	960	✓			✓		115.40
NETTLIK TOLEDO		30389173	RT UNV 10µL FL 960A/10	10 µL	10	960	✓			✓	✓	177.00
		30389177	RT UNV 10µL GL 960A/10	10 µL	10	960		✓		✓		105.80
		30389179	RT UNV 10µL GLS 960A/10	10 µL	10	960	✓	✓		✓		114.60
		30389176	RT UNV 10µL FGL 960A/10	10 µL	10	960	✓	✓		✓	✓	177.00
		30389190	RT UNV 20µL FL 960A/10	20 µL	10	960	<b>✓</b>			✓	<b>✓</b>	177.00
		30389171	RT UNV 100µL FL 960A/10	100 µL	10	960	✓			✓	✓	177.00
	Sc	30389187	RT UNV 200µL FL 960A/10	200 μL	10	960	✓			✓	<b>✓</b>	177.00
	Ē	30389188	RT UNV 200µL FLW 960A/10	200 μL	10	960	✓		✓	✓	✓	177.00
	Racked Tips	30389192	RT UNV 250µL L 960A/10	250 μL	10	960				✓		105.80
	Ra	30389196	RT UNV 250µL LW 960A/10	250 µL	10	960			<b>✓</b>	✓		105.80
		30389194	RT UNV 250µL LS 960A/10	250 μL	10	960	<b>✓</b>			✓		115.40
		30389197	RT UNV 250µL LSW 960A/10	250 μL	10	960	✓		✓	✓		142.00
		30389137	RT UNV 300µL L 768A/8	300 µL	8	768				✓		105.80
		30389138	RT UNV 300µL LS 768A/8	300 µL	8	768	<b>✓</b>			✓		115.40
		30389139	RT UNV 300µL FL 768A/8	300 µL	8	768	✓			<b>√</b>	<b>✓</b>	177.00
		30389167	RT UNV 1000µL L 768A/8	1000 µL	8	768				✓		105.80
		30389169	RT UNV 1000µL LS 768A/8	1000 μL	8	768	<b>✓</b>			<b>√</b>		115.40
William 2		30389166	RT UNV 1000µL FL 768A/8	1000 μL	8	768	✓			✓	✓	177.00
		30389284	GPS UNV 10µL GL 960A/10	10 μL	10	960		<b>√</b>		<b>√</b>		\$ 105.80
	<u>.</u>	30389286	GPS UNV 10µL GLS 960A/10	10 µL	10	960	✓	✓		✓		115.20
	Sav	30389288	GPS UNV 250µL L 960A/10	250 μL	10	960				<b>√</b>		105.80
Rack base not included	SpaceSaver	30389290	GPS UNV 250µL LS 960A/10	250 µL	10	960	✓			✓		115.20
with SpaceSaver. <u>See page 53</u> for	ß	30389280	GPS UNV 1000µL L 768A/8	1000 μL	8	768				<b>√</b>		105.80
empty racks.		30389282	GPS UNV 1000µL LS 768A/8	1000 μL	8	768	✓			✓		114.60
		17014395	RC UNV 10µL L 1000A/1	10 µL		1000				<b>✓</b>		\$ 78.60
		17014396	RC UNV 10µL GL 1000A/1	10 μL		1000		<b>√</b>		1		78.60
	Bulk	17014397	RC UNV 250µL L 1000A/1	250 μL		1000				1		78.60
	<b>m</b>	17014398	RC UNV 1000µL L 1000A/1	1000 μL		1000				1		78.60
		17015068	RC UNV 1000µL LW 1000A/1	1000 μL		1000			<b>√</b>	<b>√</b>		78.60



# **Specialty Tips**

# Unique Tips for Special Applications





#### **Extended-length tips**

For tall, narrow tube work. BioClean Ultra.

#### 102 mm long

Rainin extended-length LTS tips reach to the bottom of blocks, tall vials, flasks, 100 mm tubes, etc.

#### Narrow diameter

These tips will pass through 8 mm ID penetrable septa, such as on specimen collection kits.





#### Wide-orifice tips

Ideal for delicate samples. BioClean Ultra.

For samples with mammalian cells or high-molecular-weight DNA.
1.5 mm wide orifice minimizes shear forces on delicate samples and flow resistance for viscous samples. 200 and 1000 µL.





#### ShaftGard™ 10 µL tips

Avoid cross-contamination. BioClean Ultra.

#### Tip protects shaft, tip ejector

Shaft and tip ejector are protected from contamination, reducing the need to decontaminate or sterilize the shaft and tip ejector when working with biological, radioactive or other critical samples.

			Max V02(M.	DL	IPSCAPTON	/	_	ATURE	/		
LTS	CATALOG No.	DESCRIPTION	1 thin	PACE	20/20/20/20/20/20/20/20/20/20/20/20/20/2	STED	WIO.	0.3 MO7	FILTE	4-(F.	PRICE
	30389251	RT LTS 250µL X 768A/8	200 μL²	8	768					✓	\$ 133.40
	30389252	RT LTS 250µL SX 768A/8	200 μL²	8	768	✓				✓	143.20
	30389242	RT LTS 200µL FX 768A/8	200 μL²	8	768	✓			✓	✓	223.40
	30389222	RT LTS 1000µL X 768A/8	1000 μL	8	768					✓	133.40
	30389224	RT LTS 1000µL SX 768A/8	1000 μL	8	768	✓				✓	143.20
	30389223	RT LTS 1000µL FX 768A/8	1000 μL	8	768	✓			✓	✓	223.40
sd	30389247	RT LTS 250µL W 960A/10	200 μL²	10	960		✓				\$ 97.60
Racked Tips	30389248	RT LTS 250µL LW 960A/10	200 μL²	10	960		✓	✓			105.80
cke	30389249	RT LTS 250µL SW 960A/10	200 μL²	10	960	✓	✓				115.40
2	30389250	RT LTS 250µL LSW 960A/10	200 μL²	10	960	✓	✓	✓			142.00
	30389241	RT LTS 200µL FLW 960A/10	200 μL	10	960	✓	✓	✓	✓		177.00
	30389217	RT LTS 1000µL W 768A/8	1000 μL	8	768		✓				97.60
	30389219	RT LTS 1000µL LW 768A/8	1000 μL	8	768		✓	✓			105.80
	30389220	RT LTS 1000µL SW 768A/8	1000 μL	8	768	✓	✓				115.40
	30389221	RT LTS 1000µL LSW 768A/8	1000 μL	8	768	✓	✓	✓			142.00
	30389218	RT LTS 1000µL FLW 768A/8	1000 μL	8	768	✓	✓	✓	✓		177.00

 $<sup>^{\</sup>text{2}}$  Fits 50, 100 and 200  $\mu\text{L}$  pipettes

				DL	QTY	/	<del>/</del>	ATURES	3		
UNV	CATALOG No.	DESCRIPTION	MAY VOLUM	PACE	TIPS/CARTON	STER.	SHAFT	WIDE WIDE	0., MO7	FILTES	PRICE
	30389174	RT UNV 10µL G 960A/10	10 μL	10	960		✓				\$ 81.40
	30389178	RT UNV 10µL GS 960A/10	10 μL	10	960	✓	✓				105.80
Tips	30389175	RT UNV 10µL FG 960A/10	10 μL	10	960	✓	✓			✓	162.80
ked	30389195	RT UNV 250µL W 960A/10	250 µL	10	960			✓			89.20
Racked	30389196	RT UNV 250µL LW 960A/10	250 µL	10	960			✓	✓		105.80
	30389197	RT UNV 250µL LSW 960A/10	250 µL	10	960	✓		✓	✓		142.00
	30389188	RT UNV 200µL FLW 960A/10	250 μL	10	960	✓		✓	✓	✓	177.00
Bulk	17001127	RC UNV 250µL W 1000A/1	250 μL		1000			✓			\$ 62.60
五	17015068	RC UNV 1000µL LW 1000A/1	1000 μL		1000			✓	✓		78.60
ဟ	30389283*	GPS UNV 10µL G 960A/10	10 μL	10	960		<b>✓</b>				\$ 81.60
Refills	30389285*	GPS UNV 10µL GS 960A/10	10 μL	10	960	✓	✓				105.80
~	30389264	GP UNV 10µL FG 960A/10	10 μL	10	960	✓	✓			✓	162.80

Wide-O = Wide Orifice, Low Ret = Low Retention \*Green-Pak SpaceSaver



# The Highest Standard for Tip Purity and Performance

Rainin continuously tests and certifies its BioClean Ultra™ pipette tip products. Researchers can have absolute confidence in the integrity of their experiments and the fidelity of their data.



#### Why BioClean Ultra?

This table shows the thoroughness and completeness of Rainin BioClean Ultra quality testing, compared to other major tip brands. Rainin maintains the industry's most comprehensive testing procedures.

	HUMAN DNA	BACTERIAL DNA	DNASE
<b>Bio</b> Clean utra <sup>™</sup>	< 0.32 pg	< 1 pg	≤ 10 <sup>-7</sup> KU/µL
Brand – A	"Free of"	_	10 <sup>-5</sup> KU/μL
Brand – B	< 0.40 fg	-	-
Brand – C	< 2 pg	< 50 fg	10 <sup>-4</sup> KU
Brand — D	None detected	_	None detected
Brand — E	< 1 pg	_	"Free of"
Brand — F	"Free of"	"Free of"	"Free of"
Brand — G	-	_	< 6.25 x 10 <sup>-5</sup> U/μL
Brand — H	< 30 pg	_	<10x <sup>-7</sup> KU/µL
Brand — I	< 30 pg	_	<10 <sup>-7</sup> KU/µL
Brand – J	"Free of"	"Free of"	"Free of"



#### **Tested and Certified Purity**

- The industry's most comprehensive testing procedures.
- The only tips certified to be both proteinand protease-free.
- Virgin polypropylene completely inert, with no bioactive components.



#### Clean Manufacturing and Packaging

- Manufactured under Class 100,000 clean room conditions.
- Fully automated material handling zero human contact.
- Strict conformance to ISO 9001 eliminates all external sources of contamination.

RNASE	ENDOTOXIN	PROTEIN	PROTEASE	ATP	PCR INHIBITORS
≤ 10 <sup>-9</sup> KU/µL	≤ 0.001 EU/mL	< 2 ng/sample	≤ 500 ng/mL	< 2 x 10 <sup>-12</sup> mg/μL	None detected
10 <sup>-9</sup> KU/μL	≤ 0.05 EU/mL	-	-	-	-
< 8.6 fg	< 1 pg	-	-	< 1 fg	-
10 <sup>-9</sup> KU	< 0.001 EU/mL	<1ng/µL	_	< 5.5 x 10 <sup>-12</sup> mg	< 10 targets amplifiable
None detected	< 0.06 EU/mL	-	-	-	None detected
"Free of"	< 0.05 EU/mL	_	"Free of"	"ATP-free"	_
"Free of"	"Free of"	-	"Free of"	"ATP-free"	"Free of"
< 3.125 x 10 <sup>-9</sup> U/μL	< 0.03 EU/mL	_	_	_	_
< 10 <sup>-9</sup> KU/µL	< 0.06 EU/mL	-	-	-	_
< 10 <sup>-9</sup> KU/µL	≤ 0.06 EU/mL	_	_	$< 10^{-13}$ mg/ $\mu$ L	None detected
"Free of"	"Free of"	-	-	-	-

KU Kunitz Units EU Endotoxin Units ng nanogram (Factor: 10-9) pg picogram (Factor: 10-12) fg femtogram (Factor: 10-15)
"Non detected", "Free of" and "... Free" designations apparently indicate no observed contamination. No specifications or test methodology provided.
"-" no data available.

# Flexible, Integrated and Easy to Implement

# Pipette Management Simplified

- Assign locations, users and groups.
- Search and filter custom reports and schedules.
- Know everything about your pipettes.

- See all calibration data.
- Schedule calibrations.
- All critical pipette data located in one single database.



lanage



# **Simplify**



- The Bluetooth® wireless technology-enabled SmartStand electronically transmits pipette serial numbers, calibration and other relevant data to PipetteX.
- Rainin calibration certificates are electronically transmitted and linked to your pipette.
- Eliminate time-consuming reading of serial numbers, barcodes or stickers.

	CATALOG No.	DESCRIPTION	PRICE
_	30475913	PipetteX - Site	\$ 1,704.00
petteX	30475915	PipetteX - Unlimited	6,386.00
Pipe Soft	30478260	PipetteX - Site, Renewal	1,704.00
	30478262	PipetteX - Unlimited, Renewal	6,386.00

## **Rainin SmartStand**

# The Intelligent Pipette Stand

By reading the RFID chip that's in every Rainin XLS and XLS+ pipette, SmartStand instantly displays the compliance status for up to four pipettes — manual or electronic. Pick up a pipette and SmartStand's color LCD screen switches to a detailed view of that pipette's service and calibration status. Works seamlessly with PipetteX!



	CATALOG No.	DESCRIPTION	PRICE
ဟ	30312897	SmartStand stand with PipetteX basic software	\$ 870.00
Stands	17001255	Universal Carousel Stand for seven pipettes	173.60
	17003024	Hang-Ups, magnetic, for three pipettes	94.40
Pipette	17004992	Shelf Hang-Ups for three pipettes	94.40
•	17006638	Adapters for E4XLS with Hang-Ups, carousel	47.00

Check out our Starter Kits with SmartStand on page 15!

## Rainin SmartCheck

# Verify the Accuracy of Any Pipette

Reduce experimental errors, between-calibration compliance failures and improve pipetting technique.

Trust your pipettes! SmartCheck is a quick and easy way to verify that your pipettes are working properly. It takes very little space on the benchtop so you can use it every day.

#### Avoid expensive rework

Reduce out-of-tolerance challenges such as having to re-do expensive experiments or incur costly downstream process waste.

#### 60 seconds or less!

Pipette verification in less than a minute: SmartCheck provides a clear Pass/Fail indication following a simple process with just four gravimetric measurements.

#### Verify any brand of pipette

Verify the performance of any brand of pipette dispensing volumes between  $10-1000~\mu L$ . You can verify individual channels of multichannel pipettes too.

#### Ideal for regulated environments

SmartCheck is the only verification tool that can operate within a 21 CFR Part 11 compliant process (requires PipetteX software).




	CATALOG No.	DESCRIPTION	PRICE
eck	30564095	SmartCheck SLS1010S	\$ 1,658.00
SmartCheck	30448961	Function Test Kit	541.00
Sme	30448967	EasyScan RFID LAN set	1,535.00

# **Rainin PipetteX**

# Software to Simplify Pipette Management

PipetteX™ uses the RFID tag embedded in every Rainin XLS pipette to automatically track and manage calibration data, pipette usage and location. Our simple software solution integrates with SmartStand, SmartCheck and Rainin Service to create a powerful, easy-to-use pipette tracking and calibration management system.

- Automatic transfer of calibration documentation from service center to software
- · Controls all elements of calibration tracking and movement
- 21 CFR Part 11 compliant management system
- Integrated in-line pipette validations
- Manage all brands of pipettes

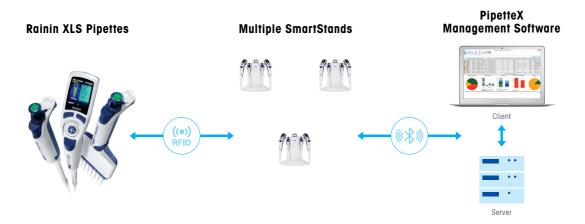

### A simple, unified system for tracking and managing pipettes (all manufacturers)

Place any Rainin XLS-brand pipette on a SmartStand to instantly see the pipette's calibration status and service details right at the bench. PipetteX connects to SmartStand via Bluetooth or USB to automatically transfer serial numbers, service and other data to the PipetteX database.

PipetteX maintains a complete history of your entire pipette inventory including calibration certificates and verification reports. This makes it easy to track the status, location, calibration and service schedule for any pipette.

That's pipette management...simplified!

mt.com/SR-PX

Prices and specifications in this publication are subject to change without notice. Rainin Terms and Conditions of Sale apply to all orders. All orders will be shipped free by UPS Ground. Alternate delivery methods are available upon request – additional charges apply. Rainin trademarks: Rainin, GPP, Good Pipetting Practice, LTS, LiteTouch, Pipet-Lite, E4, XLS, XLS+, MicroPro, BenchSmart, Liquidator, PipetteX, NanoRep, No Touch Off, SmartCheck, BioClean, BioClean Ultra, BioClean Green, Green-Pak, SpaceSaver, TerraRack, StableRak, EarthRack, ShaftGard, Rainin Classic, Pos-D, AutoRep, Disp-X, Hang-Ups. Multichannel All is a service marks of Mettler-Toledo Rainin, LLC. Eppendorf, Gilson. Integra. Thermo Fisher and Sartorius are registered trademarks of their respective owners. Bluetooth is a registered trademark of Bluetooth SiG.

Rainin MicroPro

Portable 96-Channel

**Pipetting** 

Small footprint;
 light in weight

20, 300 and
 1250 μL models

Easy to use

Affordable



New.

Streamline your 96-, 384-, and deep-well microplate work!

For more information on the new MicroPro, see page 24.

